

PC-CARD-DAS16/16

16 Analog Inputs, 16-bit A/D Resolution, Eight Digital I/O

User's Guide



PC-CARD-DAS16/16

Analog Input and Digital I/O Board

User's Guide



**MEASUREMENT
COMPUTING™**

Document Revision 1, April, 2007

© Copyright 2007, Measurement Computing Corporation

Your new Measurement Computing product comes with a fantastic extra —

Management committed to your satisfaction!

Refer to www.mccdca.com/execteam.html for the names, titles, and contact information of each key executive at Measurement Computing.

Thank you for choosing a Measurement Computing product—and congratulations! You own the finest, and you can now enjoy the protection of the most comprehensive warranties and unmatched phone tech support. It's the embodiment of our mission:

- To provide PC-based data acquisition hardware and software that will save time and save money.

Simple installations minimize the time between setting up your system and actually making measurements. We offer quick and simple access to outstanding live FREE technical support to help integrate MCC products into a DAQ system.

Lifetime warranty: Every hardware product manufactured by Measurement Computing Corporation is warranted against defects in materials or workmanship for the life of the product. Products found defective are repaired or replaced promptly.

Lifetime Harsh Environment Warranty®: We will replace any product manufactured by Measurement Computing Corporation that is damaged (even due to misuse) for only 50% of the current list price. I/O boards face some tough operating conditions, some more severe than the boards are designed to withstand. When a board becomes damaged, just return the unit with an order for its replacement at only 50% of the current list price. We don't need to profit from your misfortune. By the way, we honor this warranty for any manufacturer's board that we have a replacement for.

30 Day Money Back Guarantee: You may return any Measurement Computing Corporation product within 30 days of purchase for a full refund of the price paid for the product being returned. If you are not satisfied, or chose the wrong product by mistake, you do not have to keep it. Please call for an RMA number first. No credits or returns accepted without a copy of the original invoice. Some software products are subject to a repackaging fee.

These warranties are in lieu of all other warranties, expressed or implied, including any implied warranty of merchantability or fitness for a particular application. The remedies provided herein are the buyer's sole and exclusive remedies. Neither Measurement Computing Corporation, nor its employees shall be liable for any direct or indirect, special, incidental or consequential damage arising from the use of its products, even if Measurement Computing Corporation has been notified in advance of the possibility of such damages.

Trademark and Copyright Information

TracerDAQ, Universal Library, Harsh Environment Warranty, Measurement Computing Corporation, and the Measurement Computing logo are either trademarks or registered trademarks of Measurement Computing Corporation.

Windows, Microsoft, and Visual Studio are either trademarks or registered trademarks of Microsoft Corporation

LabVIEW is a trademark of National Instruments.

CompactFlash is a registered trademark of SanDisk Corporation.

XBee and XBee-PRO are trademarks of MaxStream, Inc.

All other trademarks are the property of their respective owners.

Information furnished by Measurement Computing Corporation is believed to be accurate and reliable. However, no responsibility is assumed by Measurement Computing Corporation neither for its use; nor for any infringements of patents or other rights of third parties, which may result from its use. No license is granted by implication or otherwise under any patent or copyrights of Measurement Computing Corporation.

All rights reserved. No part of this publication may be reproduced, stored in a retrieval system, or transmitted, in any form by any means, electronic, mechanical, by photocopying, recording, or otherwise without the prior written permission of Measurement Computing Corporation.

Notice

Measurement Computing Corporation does not authorize any Measurement Computing Corporation product for use in life support systems and/or devices without prior written consent from Measurement Computing Corporation. Life support devices/systems are devices or systems which, a) are intended for surgical implantation into the body, or b) support or sustain life and whose failure to perform can be reasonably expected to result in injury. Measurement Computing Corporation products are not designed with the components required, and are not subject to the testing required to ensure a level of reliability suitable for the treatment and diagnosis of people.

Table of Contents

Preface	
About this User's Guide	7
What you will learn from this user's guide	7
Conventions in this user's guide	7
Where to find more information	7
Chapter 1	
Introducing the PC-CARD-DAS16/16	8
Overview: PC-CARD-DAS16/16 features	8
PC-CARD-DAS16/16 block diagram	8
Software features	9
Chapter 2	
Installing the PC-CARD-DAS16/16	10
What comes with your PC-CARD-DAS16/16 shipment?	10
Hardware	10
Additional documentation	10
Optional components	10
Unpacking the PC-CARD-DAS16/16	11
Installing the software	11
Installing the PC-CARD-DAS16/16	11
If your PCMCIA card is not detected	11
Connecting the board for I/O operations	12
Connectors, cables – I/O connector	12
Pin out – I/O connector	12
Field wiring and signal termination	15
Calibrating the PC-CARD-DAS16/16	15
Chapter 3	
Programming and Developing Applications	16
Programming languages	16
Packaged applications programs	16
Register-level programming	16
Chapter 4	
Functional Details	17
Analog input resolution and range	17
Conversion speed and amplification	17
Triggering and transfer	17
A/D pacer clock	18
Chapter 5	
Specifications	19
Analog input	19
Accuracy	19
Crosstalk	20
Noise performance	20
Digital input/output	21
Interrupt	21
Counter	21

Power consumption	22
Miscellaneous	22
Environmental	22
Mechanical	22
Connector and pin out	22

About this User's Guide

What you will learn from this user's guide

This user's guide explains how to install, configure, and use the PC-CARD-DAS16/16 so that you get the most out of its analog input, digital I/O and counter features. This user's guide also refers you to related documents available on our web site, and to technical support resources.

Conventions in this user's guide

The following conventions are used in this manual to convey special information:

For more information on ...

Text presented in a box signifies additional information and helpful hints related to the subject matter you are reading.

Caution! Shaded caution statements present information to help you avoid injuring yourself and others, damaging your hardware, or losing your data.

<#:#> Angle brackets that enclose numbers separated by a colon signify a range of numbers, such as those assigned to registers, bit settings, etc.

bold text **Bold** text is used for the names of objects on the screen, such as buttons, text boxes, and check boxes. For example:
1. Insert the disk or CD and click the **OK** button.

italic text *Italic* text is used for the names of manuals and help topic titles, and to emphasize a word or phrase. For example:
The *InstaCal* installation procedure is explained in the *Quick Start Guide*.
Never touch the exposed pins or circuit connections on the board.

Where to find more information

The following electronic documents provide information relevant to the operation of the PC-CARD-DAS16/16.

- MCC's *Specifications: PC-CARD-DAS16/16* (the PDF version of the *Specifications* chapter in this guide) is available on our web site at www.mccdaq.com/pdfs/PC-CARD-DAS16_16.pdf.
- MCC's *Quick Start Guide* is available on our web site at www.mccdaq.com/PDFmanuals/DAQ-Software-Quick-Start.pdf.
- MCC's *Guide to Signal Connections* is available on our web site at www.mccdaq.com/signals/signals.pdf.
- MCC's *Universal Library User's Guide* is available on our web site at www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf.
- MCC's *Universal Library Function Reference* is available on our web site at www.mccdaq.com/PDFmanuals/sm-ul-functions.pdf.
- MCC's *Universal Library for LabVIEW™ User's Guide* is available on our web site at www.mccdaq.com/PDFmanuals/SM-UL-LabVIEW.pdf.

PC-CARD-DAS16/16 User's Guide (this document) is also available on our web site at www.mccdaq.com/PDFmanuals/PC-CARD-DAS16_16.pdf.

Introducing the PC-CARD-DAS16/16

Overview: PC-CARD-DAS16/16 features

The PC-CARD-DAS16/16 is a data acquisition and control board for IBM PC compatible computers with PCMCIA type II slots. The PC-CARD-DAS16/16 provides 16 single-ended or 8 differential analog inputs, 16-bit A/D resolution, eight digital I/O lines, and three 16-bit down counters.

The analog input range is fully programmable in one of four bipolar ranges. An on-board pacer clock, or an external pacer input, or software polling can trigger A/D conversions. Transfers are via software polling, interrupt service or REP-INSW. A 4096-word FIFO buffer provides buffering between the A/D circuit and the PCMCIA bus.

The eight digital I/O bits are available on two 4-bit ports. The digital channels are software configurable as eight inputs, eight outputs or four inputs and four outputs, and allow you to sense and control discrete events.

All signals pass through a 50-pin high-density connector. The board is completely plug-and-play, with no switches or jumpers to set.

PC-CARD-DAS16/16 block diagram

PC-CARD-DAS16/16 functions are illustrated in the block diagram shown here.

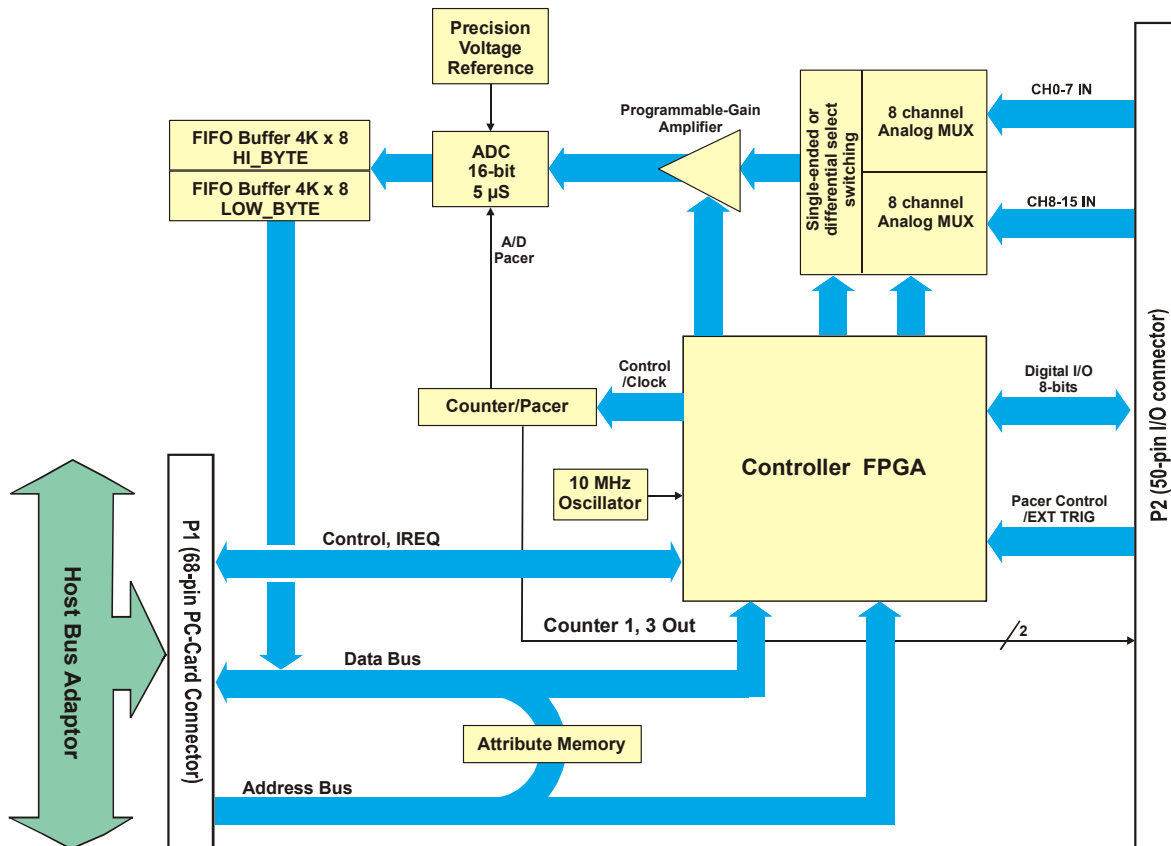


Figure 1. PC-CARD-DAS16/16 functional block diagram

Software features

For information on the features of *InstaCal* and the other software included with your PC-CARD-DAS16/16, refer to the *Quick Start Guide* that shipped with your device. The *Quick Start Guide* is also available in PDF at www.mccdaq.com/PDFmanuals/DAQ-Software-Quick-Start.pdf.

Check www.mccdaq.com/download.htm for the latest software version.

Installing the PC-CARD-DAS16/16

What comes with your PC-CARD-DAS16/16 shipment?

The following items are shipped with the PC-CARD-DAS16/16.

Hardware

- PC-CARD-DAS16/16



Additional documentation


In addition to this hardware user's guide, you should also receive the *Quick Start Guide* (available in PDF at www.mccdaq.com/PDFmanuals/DAQ-Software-Quick-Start.pdf). This booklet supplies a brief description of the software you received with your PC-CARD-DAS16/16 and information regarding installation of that software. Please read this booklet completely before installing any software or hardware.

Optional components

- Cables



CPCC-50F-39

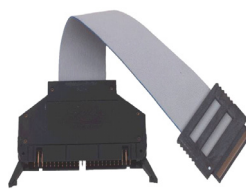

CPCC-50M-4



C50FF-x

- Signal termination and conditioning accessories

MCC provides signal conditioning and termination products for use with the PC-CARD-DAS16/16. Refer to [Field wiring and signal termination](#) on page 13 for a complete list of compatible accessory products.

Unpacking the PC-CARD-DAS16/16

As with any electronic device, you should take care while handling to avoid damage from static electricity. Before removing the PC-CARD-DAS16/16 from its packaging, ground yourself using a wrist strap or by simply touching the computer chassis or other grounded object to eliminate any stored static charge.

If any components are missing or damaged, notify Measurement Computing Corporation immediately by phone, fax, or e-mail:

- Phone: 508-946-5100 and follow the instructions for reaching Tech Support.
- Fax: 508-946-9500 to the attention of Tech Support
- Email: techsupport@mccdaq.com

Installing the software

Refer to the *Quick Start Guide* for instructions on installing the software on the *Measurement Computing Data Acquisition Software CD*. This booklet is available in PDF at www.mccdaq.com/PDFmanuals/DAQ-Software-Quick-Start.pdf.

Installing the PC-CARD-DAS16/16

The PC-CARD-DAS16/16 board is completely plug-and-play. There are no switches or jumpers to set. To install your board, follow the steps below.

Install the MCC DAQ software before you install your board

The driver needed to run your board is installed with the MCC DAQ software. Therefore, you need to install the MCC DAQ software before you install your board. Refer to the *Quick Start Guide* for instructions on installing the software.

To install your PC-Card, do the following:

- Insert the card into a free PC Card/PCMCIA type II or III slot. The key helps to insure that the cable is inserted in the correct orientation.

You do not have to turn the computer off. The system is designed for power-on installation. You should hear an insertion beep when you insert the card.



Figure 2. End view of the 50-pin PC-CARD connector showing proper orientation

Windows automatically detects, recognizes, and configures the PC-CARD. You should hear an insertion beep when you insert the card into the slot. To verify that the card is recognized, go to Control Panel\System\Device Manager and the card should now appear under "DAS Component."

If your PCMCIA card is not detected

If the card is not detected by Windows, and you are not prompted for a driver after inserting the card, check that your computer's 32-bit PCMCIA drivers are installed and enabled. Do the following:

1. From your desktop, right-click on My Computer and select **Properties**. The **System Properties** dialog opens.

2. Select the **Hardware** tab and click on the **Device Manager** button.
3. Verify that "PCMCIA adapters" is listed in the Device Manager. If you don't find this entry, or if the properties for the adapter indicate "this device is not working," you need to install or update your PCMCIA adapter drivers.
 - o If the PCMCIA adapter is not listed, use the **Add New Hardware Wizard** to install PCMCIA support.
 - o If the PCMCIA adapter is listed but not working, use the **Update Driver** option to install the appropriate drivers.

After performing the update procedure, reboot your PC and insert your card again.

Connecting the board for I/O operations

Connectors, cables – I/O connector

The table below lists the board connector, applicable cables, and compatible accessory products.

Board connector, cables, and accessory equipment

Connector type	50-pin connector
Compatible cables	<ul style="list-style-type: none"> ▪ CPCC-50F-39: 50-pin Micro connector to 50-pin female IDC, one-meter cable (39 inches). ▪ CPCC-50M-4: 50-pin Micro connector to 50-pin male IDC, 4 inch adapter cable. and <ul style="list-style-type: none"> ▪ C50FF-x: 50-pin IDC female to female cable. x = length in feet.
Compatible accessory products	CIO-MINI50 SCB-50

Pin out – I/O connector

Figure 3 shows a PC-CARD-DAS16/16 case looking into the male mini-connector. The connector is mechanically keyed to insure that the cable is inserted correctly.



Chassis Ground & Digital Ground on Connector Housing & Shield

Figure 3. 50-pin I/O mini-connector

Cabling

Measurement Computing offers two cables for connecting the PC-CARD-DAS16/16 to a screw-type terminal board or other signal conditioning interface board:

- The CPCC-50F-39 cable: 39 inches (990 mm) long; compatible with standard 50-pin screw terminal products.
- The CPCC-50M-4 cable: four-inch long adapter cable; required when using a C50FF-x series cable.

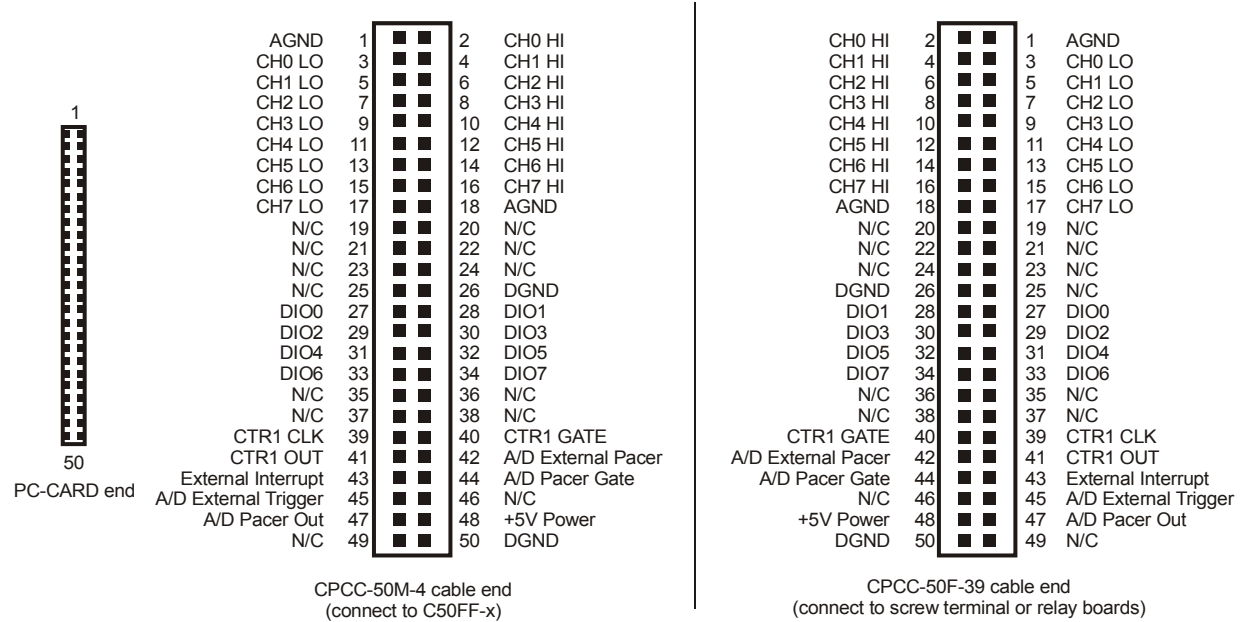


Figure 4. Differential mode cable map — PC-CARD to CPCC-50M-4 and to CPCC-50F-39

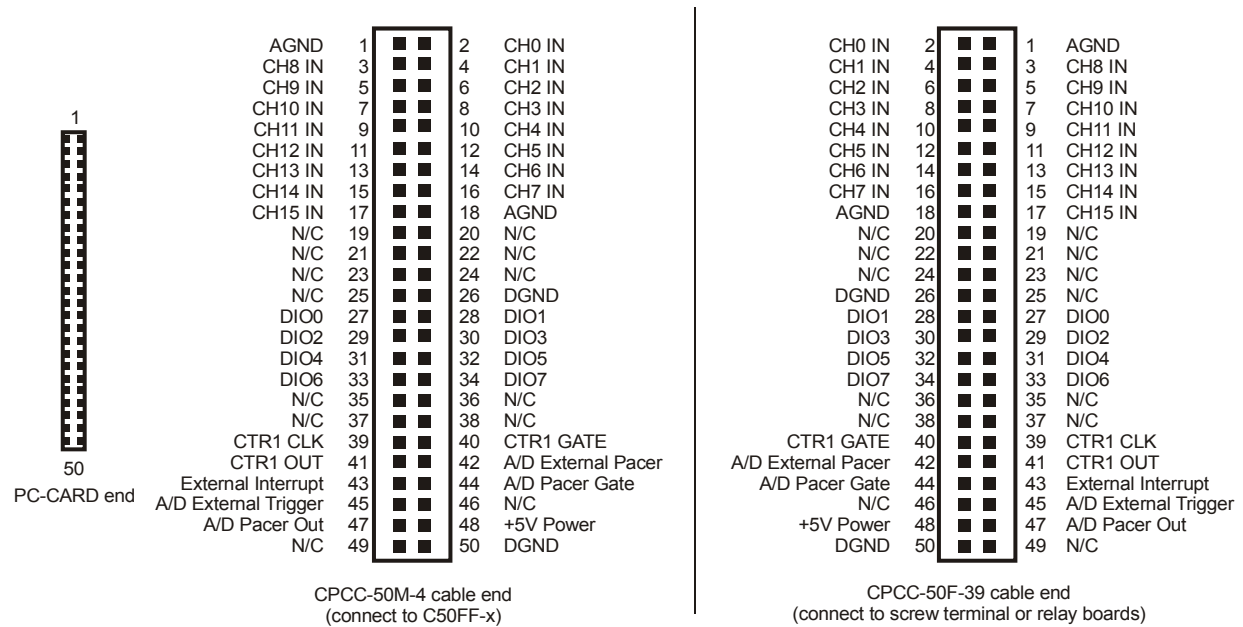


Figure 5. Single-ended mode cable map — PC-CARD to CPCC-50M-4 and to CPCC-50F-39

Note

Digital signals should not be grounded to an analog ground (AGND) pin. Use a digital ground (DGND) pin.

Caution! Do not exceed the input specifications. There are no socketed or user serviceable parts in a PC-CARD-DAS16/16. Check the specifications and input voltages *before* connecting any signals.

CPCC-50F-39

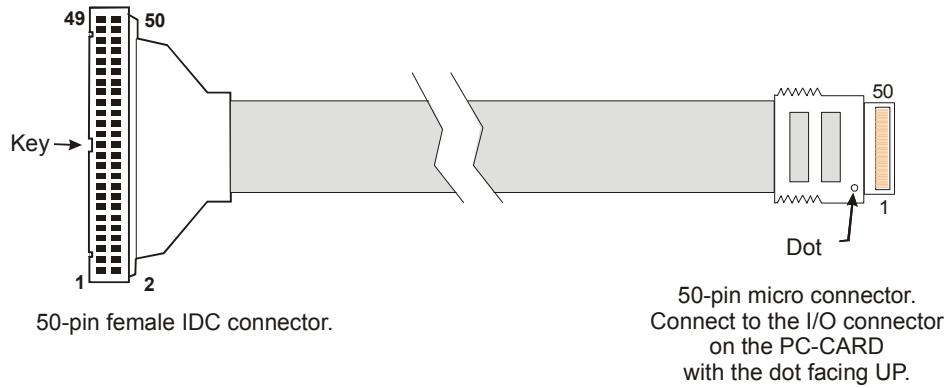


Figure 6. CPCC-50F-39 cable connections

Details on the CPCC-50F-39 cable are available on our web site at www.mccdaq.com/cbicalog/cbiproduct.asp?dept_id=105&pf_id=1379.

CPCC-50M-4

If your application requires a cable that is longer than one meter in length, use the CPCC-50M-4 four-inch cable, and connect to a C50FF-x cable.

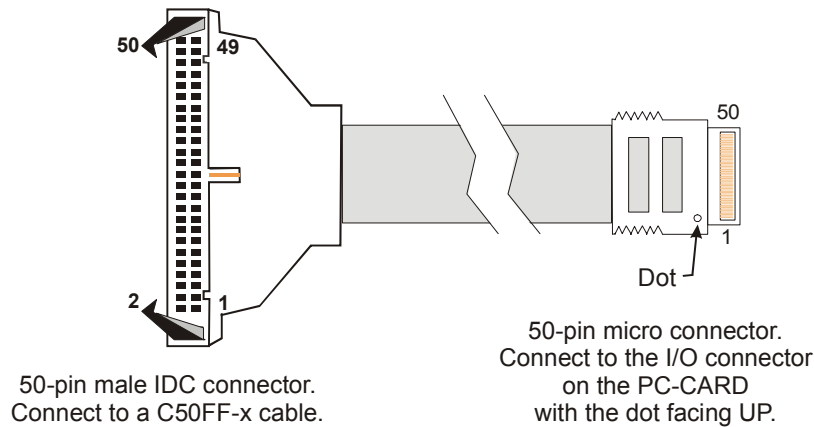


Figure 7. CPCC-50M-4 cable connections

Details on the CPCC-50M-4 cable are available on our web site at www.mccdaq.com/cbicalog/cbiproduct.asp?dept_id=96&pf_id=1380.

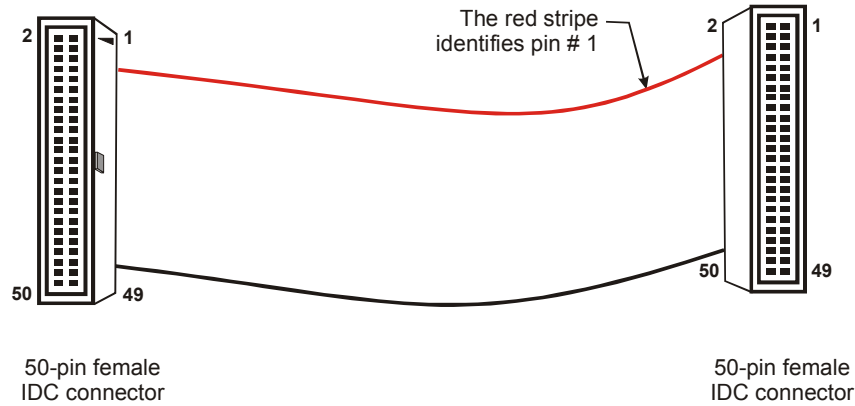
C50FF-x

Figure 8. C50FF-x cable

Details on the C50FF-x cable are available on our web site at www.mccdaq.com/cbicatalog/cbiproduct.asp?dept_id=104&pf_id=136.

Field wiring and signal termination

You can use the following cabling, screw termination, and signal conditioning products with the CPCC-50F-39 cable, or with the CPCC-50M-4 and C50FF-x cables:

- CIO-MINI50 – 50-pin screw terminal board. Details on this product are available on our web site at www.mccdaq.com/cbicatalog/cbiproduct.asp?dept_id=102&pf_id=258.
- SCB-50 – 50 conductor, shielded signal connection/screw terminal box provides two independent 50-pin connections. Details on this product are available on our web site at www.mccdaq.com/cbicatalog/cbiproduct.asp?dept_id=196&pf_id=1168.

Information on signal connections

General information regarding signal connection and configuration is available in the *Guide to Signal Connections* (available at www.mccdaq.com/signals/signals.pdf).

Calibrating the PC-CARD-DAS16/16

The PC-CARD-DAS16/16 features auto-calibration with *InstaCal*. Calibration coefficients for each range are stored on board in nonvolatile RAM. The normal calibration interval is once per year.

Programming and Developing Applications

After following the installation instructions in Chapter 2, your board should now be installed and ready for use. In general there may be no correspondence among registers for different boards. Software written at the register-level for other models does not function correctly with your board.

Programming languages

Measurement Computing's Universal Library provides access to board functions from a variety of Windows programming languages. If you are planning to write programs, or would like to run the example programs for Visual Basic® or any other language, please refer to the *Universal Library User's Guide* (available on our web site at www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf).

Packaged applications programs

Many packaged application programs now have drivers for your board. If the package you own does not have drivers for the board, please fax or e-mail the package name and the revision number from the install disks. We will research the package for you and advise how to obtain drivers.

Some application drivers are included with the Universal Library package, but not with the application package. If you have purchased an application package directly from the software vendor, you may need to purchase our Universal Library and drivers. Please contact us by phone, fax or e-mail:

- Phone: 508-946-5100 and follow the instructions for reaching Tech Support.
- Fax: 508-946-9500 to the attention of Tech Support
- Email: techsupport@mccdaq.com

Register-level programming

You should use the Universal Library or one of the packaged application programs mentioned above to control your board. Only experienced programmers should try register-level programming.

Functional Details

Analog input resolution and range

The 16-bit A/D converter provides a resolution of 1/65,536 parts of full scale. The smallest reading of full scale (1 part in 65,536) is called a Least Significant Bit (LSB). Four different bipolar ranges are controlled by software:

Analog input ranges

Bipolar	1 LSB
±10 V	0.000305 V
±5 V	0.000153 V
±2.5 V	0.000076 V
±1.25 V	0.000038 V

The input range is controlled by a programmable-gain amplifier.

Conversion speed and amplification

The A/D chip always runs at full speed. The A/D converter and sample & hold circuit captures and digitizes a signal in 5 μ s. The conversion speed of the A/D remains constant in all conditions and at all throughput rates. When you request a sample rate of say 20 kHz, the A/D converter is still converting the signal in 5 μ s. The 20 kHz rate comes from the fact that conversions are being initiated only every 50 μ s.

What factors limit conversion speed?

The first is clearly the A/D. A 5 μ s conversion speed translates to a maximum throughput of 200 kHz. The second limiting factor can be the analog front end.

The front end may consist of a multiplexer and a programmable gain amplifier. The speed at which these circuits can switch may also limit the throughput of the A/D board. That is, the rate at which it can acquire, convert and transfer a signal with full accuracy. Accuracy is the key term here. The A/D can always run at full speed, but has the front end settled and captured a true, accurate signal?

What about input range vs. speed?

The design of the analog front end is crucial to maintaining total throughput. Most A/D chips have a fixed input range, typically ± 5 V. The analog front end amplifies low level signals and adjusts unipolar signals to match the A/D converter's standard input. The PC-CARD-DAS16/16 achieves 200 kHz in all of the four ranges for single channel acquisitions.

Triggering and transfer

A trigger begins an acquisition/transfer cycle. There are three ways to trigger a PC-CARD-DAS16/16 — programmable pacer, software, or external. The trigger source selection is programmable.

The programmable pacer is the quotient of two 16-bit counters dividing a 10 MHz or 1 MHz wave derived from a 10 MHz crystal oscillator which can be used to trigger any number of paced conversions. A single conversion can be triggered by software at any time. External trigger, pacer clock and gate signals may also be used to control conversions and synchronize to external events.

Specifications

Typical for 25 °C unless otherwise specified.

Specifications in *italic text* are guaranteed by design.

Analog input

Table 1. Analog input specifications

A/D converter type	AD976A
Resolution	16 bits
Number of channels	16 single-ended / 8 differential, software selectable
Input ranges	± 10 V, ± 5 V, ± 2.5 V, ± 1.25 V, software programmable
A/D pacing (software programmable)	<ul style="list-style-type: none"> ▪ Internal counter - 82C54 ▪ External source - A/D External Pacer, software programmable for rising or falling edge ▪ Software polled
A/D trigger sources	External edge trigger (A/D External Trigger)
A/D triggering modes	Rising or falling edge trigger - software selectable
A/D gate sources	A/D External Trigger, gate high or low, software selectable
	A/D Pacer Gate, gate high
Burst Mode	Software selectable option, burst rate = 100 kHz
Data Transfer	From 4 k sample FIFO via REPINSW
	Programmed I/O
<i>A/D conversion time</i>	<i>5 μs max</i>
Calibrated throughput	200 kHz single channel, 100 kHz multiple channel. Minimum system requirement is Pentium II, 400 MHz.
Calibration	Auto-calibration, calibration factors for each range stored on board in nonvolatile RAM

Accuracy

Accuracies are listed for a 200 kHz sampling rate, single channel operation, a 60 minute warm-up, and operational temperatures within ± 2 °C of internal calibration temperature. The calibrator test source high side is tied to Channel 0 In, and the low side tied to AGND.

Table 2. Absolute Accuracy specifications

Range	Absolute accuracy
± 10.00 V	± 5.0 LSB max
± 5.000 V	± 5.0 LSB max
± 2.500 V	± 5.0 LSB max
± 1.250 V	± 5.0 LSB max

Each PC-CARD-DAS16/16 is tested at the factory to assure the board's overall error does not exceed accuracy limits described in Table 2.

Table 3. Calibrated accuracy specifications

Range	Gain Error	Offset Error	DLE (Note 1)	ILE (Note 1)
±10.00 V	±3 max	±1.5 max	-1.0, +1.75 max	±2.0 max
±5.000 V	±3 max	±1.5 max	-1.0, +1.75 max	±2.0 max
±2.500 V	±3 max	±1.5 max	-1.0, +1.75 max	±2.0 max
±1.250 V	±3 max	±1.5 max	-1.0, +1.75 max	±2.0 max

Note 1: These are the intrinsic specifications of the ADC. Software calibration may introduce a small additional amount of linearity error.

As shown in Table 3, total board error is a combination of gain, offset, differential linearity and integral linearity error. The theoretical worst-case error of the board may be calculated by summing these component errors. Worst case errors are realized only in the unlikely event that each of the component errors are at their maximum level, and causing error in the same direction.

Analog input full-scale gain drift	±0.60 LSB/°C max
Analog input zero drift	±0.15 LSB/°C max
Overall analog input drift	±0.75 LSB/°C max
Common mode range	±10 V min
CMRR @ 60 Hz	-76 dB min
Input leakage current	±20 nA max
Input impedance	10 MOhms min
Absolute maximum input voltage	+55/-40 V (Fault Protected via Input Mux)

Crosstalk

Crosstalk is defined here as the influence of one channel upon another when scanning two channels at the maximum rate. A full scale 100 Hz triangle wave is input on channel 1; channel 0 is tied to analog ground at the connector. The table below summarizes the influence of channel 1 on channel 0 with the effects of noise removed. The residue on channel zero is described in LSB's.

Table 4. Channel to channel crosstalk specifications

Condition	Crosstalk	Per channel Rate	ADC Rate
±10.00V	5LSB _{pk-pk}	50 kHz	100 kHz
±5.000V	6LSB _{pk-pk}	50 kHz	100 kHz
±2.500V	7LSB _{pk-pk}	50 kHz	100 kHz
±1.250V	10LSB _{pk-pk}	50 kHz	100 kHz

Noise performance

Table 5 summarizes the noise performance for the PC-CARD-DAS16/16. Noise distribution is determined by gathering 50 K samples at 200 kHz with inputs tied to ground at the user connector.

Table 5. Noise performance specifications

Range	% within ±2 LSBs	% within ±1 LSB	Typical LSB _{rms} *	Max LSB _{rms} *
All ranges	78%	47%	1.8	4.7

* RMS noise is defined as the peak-to-peak bin spread divided by 6.6.

Digital input/output

Table 6. DIO specifications

Digital type	FPGA
Number of I/O	8
Configuration	Two ports, four bits each. Programmable as 8 input, 8 output, or 4 input and 4 output
Input low voltage	0.8 V max
Input high voltage	2.0 V min
Output low voltage (IOL = 4 mA)	0.32 V max
Output high voltage (IOH = -4 mA)	3.86 V min
<i>Absolute maximum input voltage</i>	-0.5 V, +5.5 V
Power-up / reset state	Input mode (high impedance)

Interrupt

Table 7. Interrupt specifications

Interrupts	Programmable: Levels 2 – 15
Interrupt enable	Programmable. Default = disabled.
Interrupt sources	External (External Interrupt)
	A/D End-of-channel-scan
	A/D FIFO-not-empty
	A/D FIFO-half-full
	A/D Pacer

Counter

Table 8. Counter specifications

Counter type	82C54
Configuration	3 down counters, 16 bits each
Counter 1 - User counter	Source: Programmable external (Ctr 1 Clk) or 100 kHz internal source
	Gate: Available at connector (Ctr 1 Gate), pulled to logic high via 10K resistor (See Note 2)
	Output: Available at connector (Ctr 1 Out)
Counter 2 - ADC Pacer Lower Divider	Source: Programmable, 1MHz or 10 MHz internal source
	Gate: Available at connector (A/D Pacer Gate), pulled to logic high via 10K resistor.
	Output: Chained to Counter 3 Clock
Counter 3 - ADC Pacer Upper Divider	Source: Counter 2 Output
	Gate: Internal
	Output: Programmable as ADC Pacer clock. Available at user connector (ADC Pacer out)
<i>Clock input frequency</i>	10 MHz max
<i>High pulse width (clock input)</i>	30 ns min
<i>Low pulse width (clock input)</i>	50 ns min
<i>Gate width high</i>	50 ns min
<i>Gate width low</i>	50 ns min
<i>Input low voltage</i>	0.8 V max
<i>Input high voltage</i>	2.0 V min
<i>Output low voltage</i>	0.4 V max
<i>Output high voltage</i>	3.0 V min

Crystal oscillator frequency	10 MHz
Frequency accuracy	50 ppm

Note 2: If you are not driving the gate of User Counter 1, it is strongly recommended that it be connected to +5V (VDD).

Power consumption

Table 9. Power consumption specifications

5V quiescent	130 mA typical, 150 mA max
--------------	----------------------------

Miscellaneous

Table 10. Miscellaneous specifications

+5 Volts	Available at I/O connector (+5V Power)
	Protected by resettable fuse:
	Hold current: 350 mA max @ 20 °C still air
	Trip current: 700 mA min @ 20 °C still air
	<i>Trip and recovery time:</i> 100 mS max
	On resistance: 1.3 Ohms max

Environmental

Table 11. Environmental specifications

Operating temperature range	0 to 70 °C
Storage temperature range	-40 to 100 °C
Humidity	0 to 95% non-condensing

Mechanical

Table 12. Mechanical specifications

Card dimensions	PCMCIA type II: 85.6 mm (L) x 54.0 mm (W) x 5.0 mm (H)
-----------------	--

Connector and pin out

Table 13. Connector specifications

Connector type	50-pin connector
Compatible cables	CPCC-50F-39: 50-pin Micro connector to 50-pin female IDC, one-meter cable (39 inches).
	CPCC-50M-4: 50-pin Micro connector to 50-pin male IDC, 4 inch adapter cable. and
	C50FF-x: 50-pin IDC female to female cable. x = length in feet.
Compatible accessory products	CIO-MINI50 SCB-50

Table 14. 8-channel Differential mode pin out

Pin	Signal Name	Pin	Signal Name
1	AGND	26	DGND
2	CH0 HI	27	DIO0
3	CH0 LO	28	DIO1
4	CH1 HI	29	DIO2
5	CH1 LO	30	DIO3
6	CH2 HI	31	DIO4
7	CH2 LO	32	DIO5
8	CH3 HI	33	DIO6
9	CH3 LO	34	DIO7
10	CH4 HI	35	NC
11	CH4 LO	36	NC
12	CH5 HI	37	NC
13	CH5 LO	38	NC
14	CH6 HI	39	CTR1 CLK
15	CH6 LO	40	CTR1 GATE
16	CH7 HI	41	CTR1 OUT
17	CH7 LO	42	A/D EXTERNAL PACER
18	AGND	43	EXTERNAL INTERRUPT
19	N/C	44	A/D PACER GATE
20	N/C	45	A/D EXTERNAL TRIGGER
21	N/C	46	N/C
22	N/C	47	A/D PACER OUT
23	N/C	48	VDD +5V POWER
24	N/C	49	N/C
25	N/C	50	DGND

Table 15. 16-channel Single-ended mode pin out

Pin	Signal Name	Pin	Signal Name
1	AGND	26	DGND
2	CH0 IN	27	DIO0
3	CH8 IN	28	DIO1
4	CH1 IN	29	DIO2
5	CH9 IN	30	DIO3
6	CH2 IN	31	DIO4
7	CH10 IN	32	DIO5
8	CH3 IN	33	DIO6
9	CH11 IN	34	DIO7
10	CH4 IN	35	NC
11	CH12 IN	36	NC
12	CH5 IN	37	NC
13	CH13 IN	38	NC
14	CH6 IN	39	CTR1 CLK
15	CH14 IN	40	CTR1 GATE
16	CH7 IN	41	CTR1 OUT
17	CH15 IN	42	A/D EXTERNAL PACER
18	AGND	43	EXTERNAL INTERRUPT
19	N/C	44	A/D PACER GATE
20	N/C	45	A/D EXTERNAL TRIGGER
21	N/C	46	N/C
22	N/C	47	A/D PACER OUT
23	N/C	48	VDD +5V POWER
24	N/C	49	N/C
25	N/C	50	DGND

CE Declaration of Conformity

Manufacturer: Measurement Computing Corporation
Address: 10 Commerce Way
Suite 1008
Norton, MA 02766
USA

Category: Electrical equipment for measurement, control and laboratory use.

Measurement Computing Corporation declares under sole responsibility that the product

PC-CARD-DAS16/16

to which this declaration relates is in conformity with the relevant provisions of the following standards or other documents:

EU EMC Directive 89/336/EEC: Electromagnetic Compatibility, EN 61326 (1997) Amendment 1 (1998)

Emissions: Group 1, Class A

- EN 55011 (1990)/CISPR 11: Radiated and Conducted emissions.

Immunity: EN61326, Annex A

- IEC 1000-4-2 (1995): Electrostatic Discharge immunity, Criteria C.
- IEC 1000-4-3 (1995): Radiated Electromagnetic Field immunity Criteria B.
- IEC 1000-4-4 (1995): Electric Fast Transient Burst immunity Criteria B.
- IEC 1000-4-5 (1995): Surge immunity Criteria A.
- IEC 1000-4-6 (1996): Radio Frequency Common Mode immunity Criteria C.
- IEC 1000-4-11 (1994): Voltage Dip and Interrupt immunity Criteria A.

Tests to IEC 1000-4-8 were not required. The PC cards do not contain components that would be susceptible to magnetic fields.

Declaration of Conformity based on tests conducted by Chomerics Test Services, Woburn, MA 01801, USA in June, 2004. Test records are outlined in Chomerics Test Report #EMI3903.04.

We hereby declare that the equipment specified conforms to the above Directives and Standards.



Carl Haapaoja, Director of Quality Assurance

Measurement Computing Corporation
10 Commerce Way
Suite 1008
Norton, Massachusetts 02766
(508) 946-5100
Fax: (508) 946-9500
E-mail: info@mccdaq.com
www.mccdaq.com