## PCI-DAS6071

Analog and Digital I/O

# **User's Guide**



## **Trademark and Copyright Information**

Measurement Computing Corporation, InstaCal, Universal Library, and the Measurement Computing logo are either trademarks or registered trademarks of Measurement Computing Corporation. Refer to the Copyrights & Trademarks section on <a href="macdaq.com/legal">macdaq.com/legal</a> for more information about Measurement Computing trademarks. Other product and company names mentioned herein are trademarks or trade names of their respective companies.

© 2006 Measurement Computing Corporation. All rights reserved. No part of this publication may be reproduced, stored in a retrieval system, or transmitted, in any form by any means, electronic, mechanical, by photocopying, recording, or otherwise without the prior written permission of Measurement Computing Corporation.

#### **Notice**

Measurement Computing Corporation does not authorize any Measurement Computing Corporation product for use in life support systems and/or devices without prior written consent from Measurement Computing Corporation. Life support devices/systems are devices or systems that, a) are intended for surgical implantation into the body, or b) support or sustain life and whose failure to perform can be reasonably expected to result in injury. Measurement Computing Corporation products are not designed with the components required, and are not subject to the testing required to ensure a level of reliability suitable for the treatment and diagnosis of people.

## **Table of Contents**

Preface About this User's Guide	5
What you will learn from this user's guide	5
Conventions in this user's guide	5
Where to find more information	5
Chapter 1	
Introducing the PCI-DAS6071	6
Overview: PCI-DAS6071 features	6
Software features	6
Chapter 2	
Installing the PCI-DAS6071	
What comes with your PCI-DAS6071 shipment?	
Hardware	
Optional components	
Unpacking the board	8
Installing the software	8
Installing the hardware	8
Configuring the hardware	8
Differential input mode	
Single-ended input mode	
DAQ-Sync configuration	9
D/A external reference	
Signal connections SCSI connector	
DAQ-Sync Connector and Pinout	
Field wiring, signal termination and conditioning	14
Chapter 3	A.F.
Functional Details	
Basic architecture	
DAQ-Sync signals	
DAQ signal timing	
SCANCLK signalSTARTSCAN signal	
A/D START TRIGGER signal	
A/D STOP TRIGGER signal	
SSH signalA/D CONVERT signal	
A/D PACER GATE signal	21
A/D EXTERNAL TIME BASE signal	
A/D STOP signal	
Waveform generation timing signals	
D/A START TRIGGER signal	27
D/A CONVERT signal D/A EXTERNAL TIME BASE signal	
General-purpose counter signal timing	
CTR1 CLK signal	
CTR1 GATE signal	30
CTR1 OUT signal	

CTR2 GATE signal	31
CTR2 OUT signal	31
Chapter 4	
Calibrating the Board	
Introduction	
Calibration theory	32
Chapter 5	0.4
Specifications	
Analog input	
Settling Time	
Parametrics	
Noise Performance	
Analog output	
Analog output pacing and triggering	
Analog trigger	
Analog I/O calibration	41
Digital I/O	41
Interrupts	42
Counters	42
Configurable AUXIN<5:0>, AUXOUT<2:0> external trigger/clocks	43
DAQ-Sync inter-board triggers/clocks	44
Power consumption	44
Environmental	44
Mechanical	44
DAQ-Sync connector	45
SCSI connector	45
Declaration of Conformity	48

## **About this User's Guide**

## What you will learn from this user's guide

This user's guide describes the Measurement Computing PCI-DAS6070 data acquisition device and lists device specifications.

## Conventions in this user's guide

#### For more information

Text presented in a box signifies additional information related to the subject matter.

Caution!	Shaded caution statements present information to help you avoid injuring yourself and others, damaging your hardware, or losing your data.
<b>bold</b> text	<b>Bold</b> text is used for the names of objects on a screen, such as buttons, text boxes, and check boxes.

italic text Italic text is used for the names of manuals and help topic titles, and to emphasize a word or phrase.

## Where to find more information

Additional information about PCI-DAS6070 hardware is available on our website at <a href="www.mccdaq.com">www.mccdaq.com</a>. You can also contact Measurement Computing Corporation by phone, fax, or email with specific questions.

- Phone: 508-946-5100 and follow the instructions for reaching Tech Support
- Fax: 508-946-9500 to the attention of Tech Support
- Email: <u>techsupport</u>@mccdaq.com
- Knowledgebase: kb.mccdaq.com

If you need to program at the register level in your application, refer to the *STC Register Map for the PCI-DAS6000 Series*. This document is available at <a href="www.mccdaq.com/registermaps/RegMapSTC6000.pdf">www.mccdaq.com/registermaps/RegMapSTC6000.pdf</a>.

## **Introducing the PCI-DAS6071**

## **Overview: PCI-DAS6071 features**

This manual explains how to install and use the PCI-DAS6071 board.

The PCI-DAS6071 board provides either 32 differential or 64 single-ended analog inputs with 12-bit resolution. Input ranges are software selectable as either bipolar or unipolar.

- Bipolar input ranges:  $\pm 10V$ ,  $\pm 5V$ ,  $\pm 2.5V$ ,  $\pm 1V$ ,  $\pm 0.5V$ ,  $\pm 0.25V$ ,  $\pm 0.1V$  and  $\pm 0.05V$ .
- Unipolar input ranges: 0 to 10V, 0 to 5V, 0 to 2V, 0 to 1V, 0 to 0.5V, 0 to 0.2V and 0 to 0.1V.

The PCI-DAS6071 has eight lines of digital I/O, and provides two digital-to-analog outputs.

Each board has nine user-configurable trigger/clock/gate pins that are available at a 100-pin I/O connector. Six pins are configurable as inputs and three are configurable as outputs. Refer to Chapter 3 ("Functional Details") and Chapter 5 ("Specifications") for more information.

The PCI-DAS6071 provides triggering and synchronization capability. There are five trigger/strobes and a synchronizing clock provided on a 14-pin header. The DAQ-Sync signals use dedicated pins. Only the direction can be set. Refer to Chapter 2 ("Installing the Board") and Chapter 5 ("Specifications") for more information on these signals.

Interrupts can be generated by up to seven ADC sources and four DAC sources. Interrupt sources are listed in Chapter 5 ("Specifications").

The PCI-DAS6071 board contains an 82C54 counter chip, which consists of three 16-bit counters. Clock, gate, and output signals from two of the three counters are available on the 100-pin I/O connector. The third counter is used internally.

### Software features

For information on the features of *Insta*Cal and the other software included with your PCI-DAS6071, refer to the *Quick Start Guide* that shipped with your device. The *Quick Start Guide* is also available in PDF at www.mccdaq.com/PDFmanuals/DAQ-Software-Quick-Start.pdf.

Check www.mccdaq.com/download.htm for the latest software version.

## **Installing the PCI-DAS6071**

## What comes with your PCI-DAS6071 shipment?

The following items are shipped with the PCI-DAS6071.

### **Hardware**

PCI-DAS6071



### **Documentation**

In addition to this hardware user's guide, you should also receive the *Quick Start Guide*. This booklet provides an overview of the MCC DAQ software you received with the device, and includes information about installing the software. Please read this booklet completely before installing any software or hardware.

## **Optional components**

Cables



Signal termination and conditioning accessories
 MCC provides signal termination products for use with the PCI-DAS6071. Refer to "Field wiring, signal termination and conditioning" on page 14 for a complete list of compatible accessory products.

## Unpacking the board

As with any electronic device, you should take care while handling to avoid damage from static electricity. Before removing the PCI-DAS6071 from its packaging, ground yourself using a wrist strap or by simply touching the computer chassis or other grounded object to eliminate any stored static charge.

If any components are missing or damaged, notify Measurement Computing Corporation immediately by phone, fax, or e-mail:

- Phone: 508-946-5100 and follow the instructions for reaching Tech Support.
- Fax: 508-946-9500 to the attention of Tech Support
- Email: techsupport@mccdaq.com

For international customers, contact your local distributor. Refer to the International Distributors section on our web site at www.mccdag.com/International.

## Installing the software

Refer to the *Quick Start Guide* for instructions on installing the software on the *Measurement Computing Data Acquisition Software CD*. This booklet is available in PDF at <a href="www.mccdaq.com/PDFmanuals/DAQ-Software-Quick-Start.pdf">www.mccdaq.com/PDFmanuals/DAQ-Software-Quick-Start.pdf</a>.

## Installing the hardware

The PCI-DAS6071 board is completely plug-and-play. There are no switches or jumpers to set. Configuration is controlled by your system's BIOS. To install your board, follow the steps below.

#### Install the MCC DAQ software before you install your board

The driver needed to run your board is installed with the MCC DAQ software. Therefore, you need to install the MCC DAQ software before you install your board. Refer to the *Quick Start Guide* for instructions on installing the software.

- 1. Turn your computer off, open it up, and insert your board into any available PCI slot.
- 2. Close your computer and turn it on.
  - If you are using an operating system with support for plug-and-play, a dialog box notifies you that new hardware has been detected. If the information file for this board is not already loaded onto your PC, you will be prompted for the disk containing this file. The MCC DAQ software contains this file. If required, insert the MCC DAQ CD and click **OK**.
- 3. To test your installation and configure your board, run the *InstaCal* utility installed in the previous section. Refer to the *Quick Start Guide* that came with your board for information on how to initially set up and load InstaCal.

Allow your computer to warm up for at least 15 minutes before acquiring data with this board. The high speed components used on the board generates heat, and it takes this amount of time for a board to reach steady state if it has been powered off for a significant amount of time.

## Configuring the hardware

All hardware configuration options on the PCI-DAS6071 are software controlled. You can select some of the configuration options using *Insta*Cal, such as the analog input configuration (64 single-ended or 32 differential channels), the edge used for triggering when using an external pacer, and the source for the two independent counters. Once selected, any program that uses the Universal Library will initialize the hardware according to these selections.

Following is an overview of the available hardware configuration options for this board. There is additional general information regarding analog signal connection and configuration in the *Guide to Signal Connections* (available on our web site at <a href="https://www.mccdaq.com/signals/signals.pdf">www.mccdaq.com/signals/signals.pdf</a>).

## Differential input mode

When all channels are configured for differential input mode, 32 analog input channels are available. In this mode, the input signal is measured with respect to the low input. The input signal is delivered through three wires:

- The wire carrying the signal to be measured connects to CH# IN HI.
- The wire carrying the reference signal connects to CH# IN LO.
- The third wire is connected to LLGND.

Differential input mode is the preferred configuration for applications in noisy environments, or when the signal source is referenced to a potential other than PC ground.

## Single-ended input mode

When all channels are configured for single-ended input mode, 64 analog input channels are available. In this mode, the input signal is referenced to the board's signal ground (LLGND). The input signal is delivered through two wires:

- The wire carrying the signal to be measured connects to CH# IN HI.
- The other wire is connected to LLGND.

## Non-referenced single-ended input mode

This mode is a compromise between differential and single-ended modes. It offers some of the advantages of each mode. Using non-referenced single-ended mode, you can still get noise rejection, but not the limitation in the number of channels resulting from a fully differential configuration. The possible downside is that the external reference input must be the same for every channel. It is equivalent to configuring the inputs for differential mode and then tying all of the low inputs together and using that node as the reference input.

When configured for non-referenced single-ended input mode, 64 analog input channels are available. In this mode, each input signal is not referenced to the board's ground, but to a common reference signal (AISENSE). The input signal is delivered through three wires:

- The wire carrying the signal to measure connects to CH# IN HI.
- The wire carrying the reference signal connects to AISENSE.
- The third wire is connected to LLGND.

This mode is useful when the application calls for differential input mode, but the limitation on channel count prevents it.

## **DAQ-Sync configuration**

You can interconnect multiple boards in the PCI-DAS6000 series to synchronize data acquisition or data output. To do this, order and install a CDS-14-*x* cable at the DAQ-Sync connectors (P2) between the boards to be synchronized.

The "x" in the CDS-14-x part number specifies the number of connectors available on the cable, and therefore, the number of boards you can interconnect. Using a CDS-14-2, you can connect two PCI-DAS6000 series boards together for I/O synchronization. Using a CDS-14-3, you can synchronize three boards, and so on. You can connect up to five PCI-DAS6000 series boards. A CDS-14-3 cable is shown in Figure 3 on page 14.

By default, all DAQ-Sync connectors are configured as inputs (slave mode). In order to be useful, one board must be set through software to serve as the master, and the signal sources of the slave boards must be defined.

#### D/A external reference

This product is equivalent in the specification and performance to the National Instruments (NI) PCI-6071E. The Analog Output External Reference Bandwidth is specified at 1 MHz. This signal path is indeed 1 MHz; however the performance of the D/A converter itself in multiplying mode varies with frequency and code. To the best of our knowledge, this behavior is not fully explained in documentation presently available from NI.

If you are using the external reference feature of the analog outputs, please refer to the data sheet for the Analog Devices AD7945. This document summarizes the performance of the D/A converter in multiplying mode.

When using a sinusoidal input signal for the Analog Output External Reference, the -3dB multiplying bandwidth of both VDAC0 and VDAC1 is approximately 250 kHz. It should also be noted that the AC feedthrough of both VDAC0 and VDAC1 increases with input frequency and may be as much as -45 dB for a 1 MHz sinusoidal External Reference input.

## **Signal connections**

## **SCSI** connector

The table below lists the board connectors, applicable cables and compatible accessory boards.

#### Board connectors, cables, accessory equipment

Connector type	Shielded SCSI 100 D-Type			
Compatible cables	C100HD50-x, unshielded ribbon cable. $x = 3$ or 6 feet			
	(see Figure 1 on page 13)			
	C100MMS-x, shielded round cable. $x = 1, 2, \text{ or } 3 \text{ meters}$			
	(see Figure 2)			
Compatible accessory products	ISO-RACK16/P			
(with the C100HD50-x cable)	ISO-DA02/P			
	BNC-16SE			
	BNC-16DI			
	CIO-MINI50			
	CIO-TERM100			
	SCB-50			
Compatible accessory products	SCB-100			
(with the C100MMS-x cable)				

32-channel differential mode pinout

Signal Name	Pin		Pin	Signal Name
GND	100	۲.,	50	GND
CTR2 OUT	99		49	AUXIN5 / A/D PACER GATE
CTR2 GATE	98		48	AUXIN4 / D/A START TRIGGER
CTR2 CLK	97			AUXIN3 / D/A UPDATE
GND	96		46	AUXIN2 / A/D STOP TRIGGER
CTR1 OUT	95		45	AUXIN1 / A/D START TRIGGER
CTR1 GATE	94		44	D/A EXTREF
CTR1 CLK	93	• •	43	AUXINO / A/D CONVERT / ATRIG
DIO7	92	• •	42	AUXOUT2 / SCANCLK
DIO6	91	• •	41	AUXOUT1 / A/D PACER OUT
DIO5	90		40	AUXOUT0 / D/A PACER OUT
DIO4	89		39	PC +5 V
DIO3	88	• •	38	D/A OUT1
DIO2	87	• •	37	D/A GND
DIO1	86	• •	36	D/A OUT0
DIO0	85	• •	35	AISENSE
CH31 IN LO	84	• •	34	CH15 IN LO
CH31 IN HI	83	• •	33	CH15 IN HI
CH30 IN LO	82	• •	32	CH14 IN LO
CH30 IN HI	81	• •	31	CH14 IN HI
CH29 IN LO	80	• •	30	CH13 IN LO
CH29 IN HI	79	• •	29	CH13 IN HI
CH28 IN LO	78	• •	28	CH12 IN LO
CH28 IN HI	77	• •	27	CH12 IN HI
CH27 IN LO	76	• •	26	CH11 IN LO
CH27 IN HI	75	• •	25	CH11 IN HI
CH26 IN LO	74	• •	24	CH10 IN LO
CH26 IN HI	73	• •	23	CH10 IN HI
CH25 IN LO	72	• •	22	CH9 IN LO
CH25 IN HI	71	• •	21	CH9 IN HI
CH24 IN LO	70	• •	20	CH8 IN LO
CH24 IN HI	69	• •	19	CH8 IN HI
LLGND	68	• •	18	LLGND
CH23 IN LO	67	• •	17	CH7 IN LO
CH23 IN HI	66	• •	16	CH7 IN HI
CH22 IN LO	65	• •	15	CH6 IN LO
CH22 IN HI	64	• •	14	CH6 IN HI
CH21 IN LO	63	• •	13	CH5 IN LO
CH21 IN HI	62	• •	12	CH5 IN HI
CH20 IN LO	61	• •	11	CH4 IN LO
CH20 IN HI	60	• •	10	CH4 IN HI
CH19 IN LO	59	• •	9	CH3 IN LO
CH19 IN HI	58	• •	8	CH3 IN HI
CH18 IN LO	57	• •	7	CH2 IN LO
CH18 IN HI	56	• •	6	CH2 IN HI
CH17 IN LO	55	••	5	CH1 IN LO
CH17 IN HI	54	• •	4	CH1 IN HI
CH16 IN LO	53	• •	3	CH0 IN LO
CH16 IN HI	52	••	2	CH0 IN HI
LLGND	51	$\overline{}$	1	LLGND
PC	l slot ↓			

64-channel single-ended mode pinout

Signal Name	Pin		Pin	Signal Name
GND	100		50	GND
CTR2 OUT	99	••	49	AUXIN5 / A/D PACER GATE
CTR2 GATE	98	••	48	AUXINA / D/A START TRIGGER
CTR2 GATE	97		47	AUXIN3 / D/A UPDATE
GND	96		46	AUXIN2 / A/D STOP TRIGGER
CTR1 OUT	95		45	AUXIN1 / A/D START TRIGGER
CTR1 GATE	94		44	D/A EXTREF
CTR1 CLK	93		43	AUXINO / A/D CONVERT / ATRIG
DIO7	92		42	AUXOUT2 / SCANCLK
DIO6	91		41	AUXOUT1 / A/D PACER OUT
DIO5	90	• •	40	AUXOUTO / D/A PACER OUT
DIO4	89	• •	39	PC +5 V
DIO3	88	• •	38	D/A OUT1
DIO2	87	• •	37	D/A GND
DIO1	86	• •	36	D/A OUT 0
DIO0	85	• •	35	AISENSE
CH63 IN	84	• •	34	CH47 IN
CH31 IN	83	• •	33	CH15 IN
CH62 IN	82	• •	32	CH46 IN
CH30 IN	81	• •	31	CH14 IN
CH61 IN	80	• •	30	CH45 IN
CH29 IN	79	• •	29	CH13 IN
CH60 IN	78	• •	28	CH44 IN
CH28 IN	77	• •	27	CH12 IN
CH59 IN	76	• •	26	CH43 IN
CH27 IN	75	• •	25	CH11 IN
CH58 IN	74	• •	24	CH42 IN
CH26 IN	73	• •	23	CH10 IN
CH57 IN	72	• •	22	CH41 IN
CH25 IN	71	• •	21	CH9 IN
CH56 IN	70	• •	20	CH40 IN
CH24 IN	69	• •	19	CH8 IN
LLGND	68	• •	18	LLGND
CH55 IN	67	• •	17	CH39 IN
CH23 IN	66	• •	16	CH7 IN CH38 IN
CH54 IN CH22 IN	65 64	••	15 14	CH38 IN
CH22 IN	63	••	13	CH37 IN
CH33 IN	62		12	CH57 IN
CH52 IN	61		11	CH36 IN
CH20 IN	60		10	CH4 IN
CH51 IN	59		9	CH35 IN
CH19 IN	58		8	CH3 IN
CH50 IN	57	••	7	CH34 IN
CH18 IN	56		6	CH2 IN
CH49 IN	55		5	CH33 IN
CH17 IN	54		4	CH1 IN
CH48 IN	53	• •	3	CH32 IN
CH16 IN	52	• •	2	CH0 IN
LLGND	51	••	1	LLGND
PC	l slot ↓			
		/		

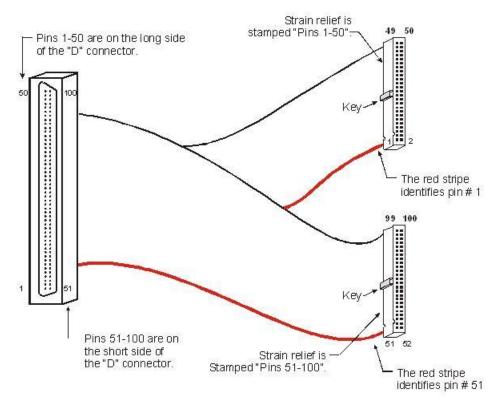


Figure 1. C100HD50-x cable

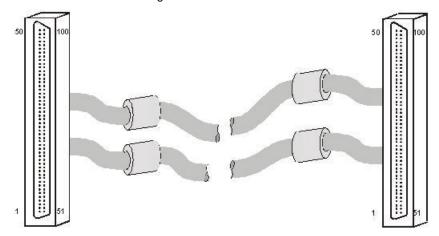


Figure 2. C100MMS-x cable

Details on these cables are available on our web site at <a href="www.mccdaq.com/products/accessories.aspx">www.mccdaq.com/products/accessories.aspx</a>.

## **DAQ-Sync Connector and Pinout**

DAQ-Sync connector and cable types

Connector type	14-pin right-angle 100 mil box header
Compatible cable	MCC p/n: CDS-14-x, 14 pin ribbon cable for board-to board DAQ-Sync connection;
_	x = number of boards (Figure 3 shows a CDS-14-3 cable)

Signal Name	Pin	 Pin	Signal Name
DS A/D START TRIGGER	1	 2	GND
DS A/D STOP TRIGGER	3	 4	GND
DS A/D CONVERT	5 .	 6	GND
DS D/A UPDATE	7	 8	GND
DS D/A START TRIGGER	9	 10	GND
RESERVED	11	 12	GND
SYNC CLK	13	 14	GND

#### DAQ-sync connector pinout (view from top)

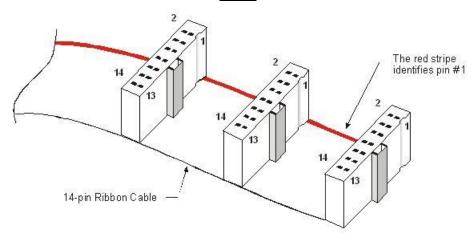


Figure 3. CDS-14-3 cable

Details on the CDS-14-x cable are available on our web site at www.mccdaq.com/cbicatalog/cbiproduct.asp?dept\_id=104&pf\_id=1528.

## Field wiring, signal termination and conditioning

#### Screw terminal boards and BNC adapters

Use with the C100HD50-x cable:

- CIO-MINI50: 50-pin universal screw terminal accessory.
- CIO-TERM100: 16x4 screw terminal.
- SCB-50: 50 conductor, shielded signal connection/screw terminal box provides two independent 50-pin connections.

Use with the C100MMS-x cable:

 SCB-100: 100 conductor, shielded signal connection/screw terminal box provides two independent 50-pin connections.

#### **BNC** connector boxes:

- BNC-16SE: 16-channel single-ended BNC connector box.
- BNC-16DI: Eight-channel differential BNC connector box.

Details on these products are available on our web site at <a href="www.mccdaq.com/products/screw">www.mccdaq.com/products/screw</a> terminal bnc.aspx.

#### ISO-5B module racks

Use with the C100HD50-x cable:

- ISO-RACK16/P: 16-channel isolation module mounting rack.
- ISO-DA02/P: Two-channel, 5B module rack.

Details on these products are available on our web site at <a href="www.mccdaq.com/products/signal\_conditioning.aspx">www.mccdaq.com/products/signal\_conditioning.aspx</a>.

## **Functional Details**

## **Basic architecture**

Figure 4 on page 17 shows a simplified block diagram of the PCI-DAS6071. This board provides all of the functional elements shown in the figure.

The System Timing and Control (STC) is the logical center for all DAQ, DIO, and DAC (if applicable) operations. It communicates over two major busses: a local bus and a memory bus.

The local bus carries digital I/O data and software commands from the PCI Bus Master. There are two Direct Memory Access (DMA) channels provided for data transfers to the PC.

Primarily, the memory bus carries A/D and D/A related data and commands. There are three buffer memories provided on the memory bus:

- The queue buffer (8K configuration memory) stores programmed channel numbers, gains, and offsets.
- The ADC buffer (8K FIFO [First In, First Out]) temporarily stores scanned and converted analog inputs.
- The *DAC 16K buffer* stores data to be output as analog waveforms.

## **Auxiliary input & output interface**

The board's 100-pin I/O connector provides six software-selectable inputs, and three software-selectable outputs. The signals are user-configurable clocks, triggers and gates.

Refer to "DAQ signal timing" on page 18 for information about these signals and their timing requirements.

The following table lists the available signals and the default signals that can be configured on the nine pins.

#### Auxiliary I/O Signals

I/O Type	Signal Name	Function		
AUXIN<5:0> sources	A/D CONVERT	External ADC Convert Strobe (default)		
(software selectable)	A/D TIMEBASE IN	External ADC Pacer Time Base		
	A/D START TRIGGER	ADC Start Trigger (default)		
	A/D STOP TRIGGER	ADC Stop Trigger (default)		
	A/D PACER GATE	External ADC Gate (default)		
	D/A START TRIGGER	DAC Trigger/Gate (default)		
	D/A UPDATE	DAC Update Strobe (default)		
	D/A TIMEBASE IN	External DAC Pacer Time Base		
AUXOUT<2:0> sources	STARTSCAN	A pulse indicating the start of conversion.		
(software selectable)	SSH	An active signal that terminates at the start of the last conversion in a scan.		
	A/D STOP	Indicates the end of a scan		
	A/D CONVERT	ADC convert pulse (default)		
	SCANCLK	Delayed version of ADC convert (default)		
	CTR1 CLK	CTR1 clock source		
	D/A UPDATE	D/A update pulse (default)		
	CTR2 CLK	CTR2 clock source		
	A/D START TRIGGER	ADC Start Trigger Out		
	A/D STOP TRIGGER	ADC Stop Trigger Out		
	A/D PACER GATE	External ADC gate		
	D/A START TRIGGER	DAC Start Trigger Out		
Default selections	AUXIN0	A/D CONVERT		
summary	AUXIN1	A/D START TRIGGER		
	AUXIN2	A/D STOP TRIGGER		
	AUXIN3	D/A UPDATE		
	AUXIN4	D/A START TRIGGER		
	AUXIN5	A/D PACER GATE		
	AUXOUT0	D/A UPDATE		
	AUXOUT1	A/D CONVERT		
	AUXOUT2	SCANCLK		

## **DAQ-Sync signals**

The DAQ-Sync hardware provides the capability of triggering or clocking up to four slave boards from a master board to synchronize data input and/or output.

The PCI-DAS6071 board provides the capability of inter-board synchronization between boards in the PCI-DAS6000 family. There are five trigger/strobes and a synchronizing clock provided on a 14-pin header. The following signals are available:

- DS A/D START TRIGGER
- DS A/D STOP TRIGGER
- DS A/D CONVERT
- DS D/A UPDATE
- DS D/A START TRIGGER
- SYNC CLK

Except for the SYNC CLK signal, the DAQ-Sync timing and control signals are a subset of the AUXIO signals available at the 100-pin I/O connector. These versions of the signals are used for board-to-board synchronization and have the same timing specifications as their I/O connector counterparts. Refer to "DAQ signal timing" on page 18 for explanations of signals and timing.

Use the SYNC CLK signal to determine the master/slave configuration of a DAQ-Sync-enabled system. Each system can have one master and up to three slaves. SYNC CLK is the 40 MHz time-base used to derive all board timing and control. The master provides this clock to the slave boards so that all boards in the DAQ-sync-enabled system are timed from the same clock.

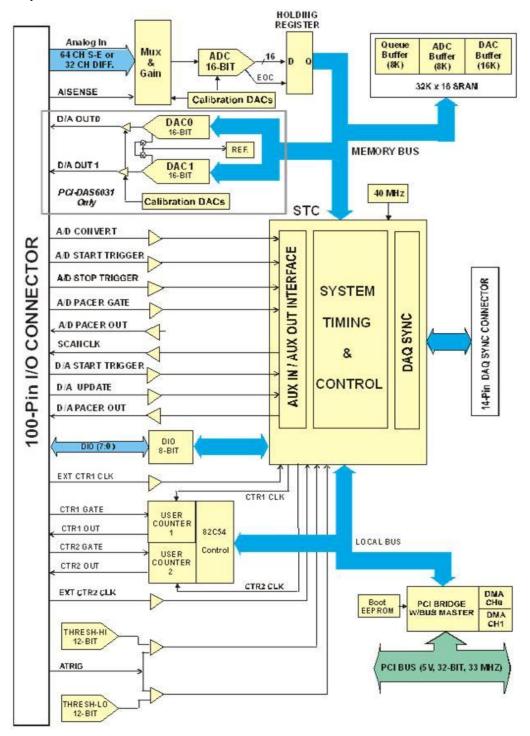


Figure 4. Functional block diagram

## **DAQ** signal timing

The DAQ timing signals are:

- SCANCLK
- A/D START TRIGGER
- A/D STOP TRIGGER
- STARTSCAN
- SSH
- A/D CONVERT
- A/D PACER GATE
- A/D EXTERNAL TIME BASE
- A/D STOP
- ATRIG

## **SCANCLK** signal

SCANCLK is an output signal that may be used for switching external multiplexers. It is a 400 ns wide pulse that follows the CONVERT signal after a 50 ns delay. This is adequate time for the analog input signal to be acquired so that the next signal may be switched in. The polarity of the SCANCLK signal is programmable. The default output pin for the SCANCLK signal is AUXOUT2, but any of the AUXOUT pins may be programmed as a SCANCLK output.

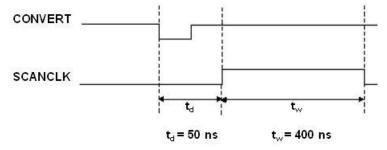


Figure 5. SCANCLK signal timing

### STARTSCAN signal

The STARTSCAN output signal indicates when a scan of channels has been initiated. You can program this signal to be available at any of the AUXOUT pins. The STARTSCAN output signal is a 50 ns wide pulse the leading edge of which indicates the start of a channel scan. Figure 6 shows the timing for the STARTSCAN signal.

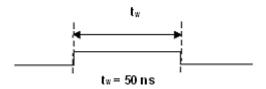


Figure 6. STARTSCAN start of scan timing

## A/D START TRIGGER signal

Use the A/D START TRIGGER signal for conventional triggering (when you only need to acquire data after a trigger event). Figure 7 shows the A/D START TRIGGER signal timing for a conventionally triggered acquisition.

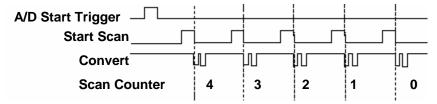


Figure 7. Data acquisition example for conventional triggering

The A/D START TRIGGER source is programmable and may be set to any of the AUXIN inputs or to the DAQ-Sync DS A/D START TRIGGER input. The polarity of this signal is also programmable to trigger acquisitions on either the positive or negative edge.

The A/D START TRIGGER signal is also available as an output and can be programmed to appear at any of the AUXOUT outputs. Refer to Figure 8 and Figure 9 for A/D START TRIGGER input and output timing requirements.

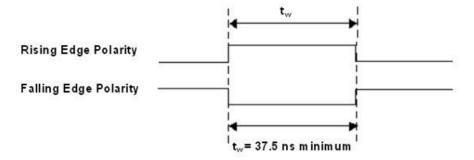


Figure 8. A/D START TRIGGER input signal timing

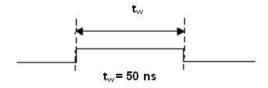


Figure 9. A/D START TRIGGER output signal timing

The A/D START TRIGGER signal is also used to initiate pre-triggered DAQ operations (when you need to acquire data just before a trigger event). In most pre-triggered applications, the A/D START TRIGGER signal is generated by a software trigger. The use of A/D START TRIGGER and A/D STOP TRIGGER in pre-triggered DAQ applications is explained next.

## A/D STOP TRIGGER signal

Pre-triggered data acquisition continually acquires data into a circular buffer until a specified number of samples have been collected after the trigger event. Figure 10 illustrates a typical pre-triggered DAQ sequence.

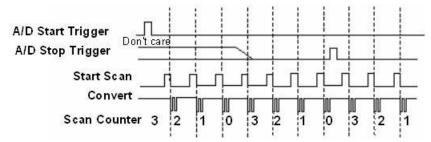


Figure 10. Pre-triggered data acquisition example

The A/D STOP TRIGGER signal signifies when the circular buffer should stop and when the specified number of post trigger samples should be acquired. It is available as an output and an input. By default, it is available at AUXIN2 as an input but may be programmed for access at any of the AUXIN pins or the DAQ-Sync "DS A/D STOP TRIGGER" input. It may be programmed for access at any of the AUXOUT pins as an output.

When using the A/D STOP TRIGGER signal as an input, the polarity may be configured for either rising or falling edge. The selected edge of the A/D STOP TRIGGER signal initiates the post-triggered phase of a pre-triggered acquisition sequence.

As an output, the A/D STOP TRIGGER signal indicates the event separating the pre-trigger data from the post-trigger data. The output is an active high pulse with a pulse width of 50 ns. Figure 11 and Figure 12 show the input and output timing requirements for the A/D STOP TRIGGER signal.

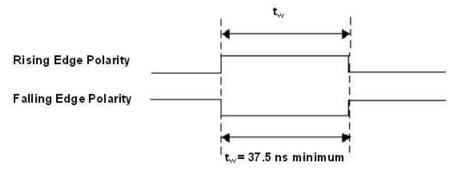


Figure 11. A/D STOP TRIGGER input signal timing

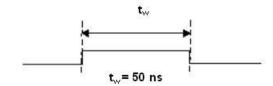


Figure 12. A/D STOP TRIGGER output signal timing

## SSH signal

The SSH signal can be used as a control signal for external sample/hold circuits. The SSH signal is a programmable polarity pulse that is asserted throughout a channel scan. The state of this signal changes after the start of the last conversion in the scan. The SSH signal may be routed via software selection to any of the AUXOUT pins. Figure 13 shows the timing for the SSH signal.

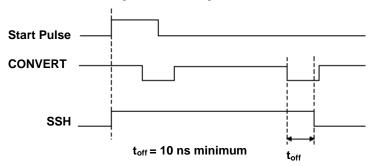


Figure 13. SSH signal timing

## A/D CONVERT signal

The A/D CONVERT signal indicates the start of an A/D conversion. It is available through software selection as an input to any of the AUXIN pins (defaulting to AUXIN0) or the DAQ-Sync DS A/D CONVERT input and as an output to any of the AUXOUT pins.

When used as an input, the polarity is software selectable. The A/D CONVERT signal starts an acquisition on the selected edge. The selected edge (either rising of falling) of the convert pulses must be separated by a minimum of 0.8 µs to remain within the 1.25 MS/s conversion rate specification.

Refer to Figure 7 and Figure 10 for the relationship of A/D CONVERT to the DAQ sequence. Figure 14 and Figure 15 show the input and output pulse width requirements for the A/D CONVERT signal.

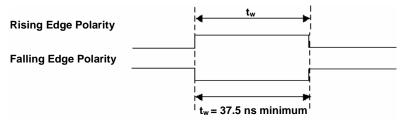


Figure 14. A/D CONVERT signal input timing requirement

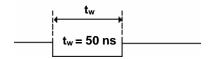


Figure 15. A/D CONVERT signal output timing requirement

The A/D CONVERT signal is generated by the on-board pacer circuit unless the external clock option is in use. This signal may be gated by hardware (A/D PACER GATE) or software.

## A/D PACER GATE signal

The A/D PACER GATE signal is used to disable scans temporarily. This signal may be programmed for input at any of the AUXIN pins.

If the A/D PACER GATE signal is active, no scans can occur. If the A/D PACER GATE signal becomes active during a scan in progress, the current scan is completed and scans are then held off until the gate is de-asserted.

## A/D EXTERNAL TIME BASE signal

The A/D EXTERNAL TIME BASE signal can serve as the source for the on-board pacer circuit rather than using the 40 MHz internal time base. Any AUXIN pin can be set programmatically as the source for this signal. The polarity is programmable.

The maximum frequency for the A/D EXTERNAL TIME BASE signal is 20 MHz. The minimum pulse width is 23 ns high or low. There is no minimum frequency specification. Figure 16 shows the timing specifications for the A/D EXTERNAL TIME BASE signal.

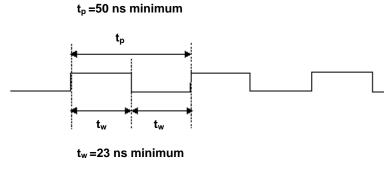


Figure 16. A/D EXTERNAL TIME BASE signal timing

## A/D STOP signal

The A/D STOP signal indicates a completed acquisition sequence. You can program this signal to be available at any of the AUXOUT pins. The A/D STOP output signal is a 50 ns wide pulse whose leading edge indicates a DAQ done condition. Figure 17 shows the timing for the A/D STOP signal.

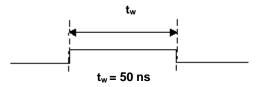


Figure 17. A/D STOP signal timing

### ATRIG signal

In addition to standard digital trigger features, the PCI-DAS6071 also provides analog triggering capability. When using the analog trigger, acquisitions may be started and controlled via an analog signal. There are four trigger/gate modes available using the analog trigger feature:

- Trigger positive or negative slope
- Gate above reference or below reference
- Hysteresis positive or negative hysteresis
- Window inside or outside window

The Trigger mode is used to start an acquisition sequence. The remaining modes provide gating functions during an acquisition sequence which start and stop the acquisition based on the gate condition.


There are two possible inputs for the analog trigger source (see Figure 18):

■ The first is the AUXINO/ATRIG pin on the 100-pin I/O connector. This is a software selectable dual-purpose pin that supports either digital or analog trigger inputs. The source selection defaults to analog trigger on power-up and may be modified at any time using *InstaCal*. The input range on the ATRIG pin is always ±10V. 8-bit DACs are used to set the HI and LO levels for the threshold(s). The threshold resolution in this mode is 78mV per step.

**Caution!** Remove all analog inputs before configuring this pin as a digital input. Any voltage levels above ±15V in this configuration may cause damage to the product!

The second possible analog trigger source is the post-gain version of any one of the 64 analog inputs. In this mode, the voltage present on the first channel in the scan may be used to initiate the acquisition sequence.

Since the input to the analog trigger circuit has been scaled by the selected range, the effective resolution of the thresholds is equal to the A/D's full-scale-range ( $\pm 2.5$ V) divided by 4096. For example, the  $\pm 2.5$ V range allows for 5V/4096, or 1.2 mV of threshold resolution.

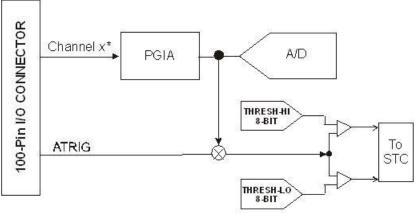

<sup>\*</sup> Only applies to the first channel in the scan

Figure 18. ATRIG circuit

The next section includes a detailed description of each mode of operation. In each case, a  $\pm 2V$  triangle waveform is used as the ATRIG input source. The THRESH\_HI is set to 1.0V and the THRESH\_LO signal is set to -1.0V.

In the following analog trigger signal diagrams, the bold portion of the waveform indicates the data acquired for the given ATRIG mode.

### **Trigger Above**

The acquisition will begin when the ATRIG signal first goes above the THRESH\_HI. This mode is non-retriggerable.

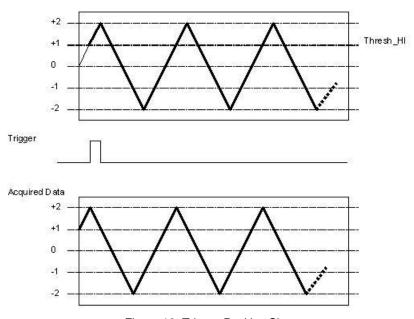


Figure 19. Trigger Positive Slope

### **Trigger Below**

The acquisition will begin when ATRIG signal first goes below the THRESH\_LO level. This mode is non-retriggerable.

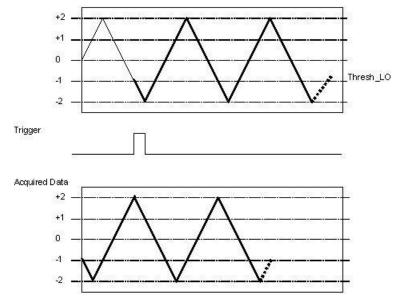


Figure 20. Trigger Negative Slope

## **Gate Above**

Data acquisition is enabled whenever ATRIG goes above the THRESH\_HI level. Acquisition is suspended whenever the ATRIG signal goes below the THRESH\_HI level. This is a level-sensitive gating mode.

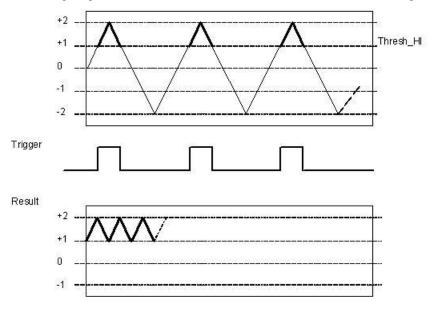


Figure 21. Gate Above

### **Gate Below**

Data acquisition is enabled whenever ATRIG goes below the THRESH\_LO level. Acquisition is suspended whenever the ATRIG signal goes above the THRESH\_LO level. This is a level-sensitive gating mode.

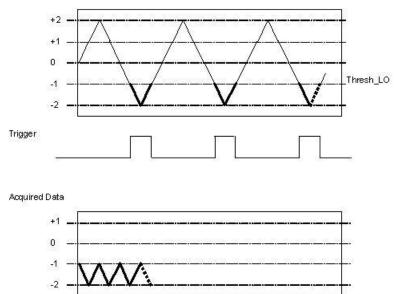


Figure 22. Gate Below

### **Gate Negative Hysteresis**

Data acquisition is enabled whenever ATRIG goes above the THRESH\_HI level. Acquisition is suspended whenever the ATRIG signal goes below the THRESH\_LO level. The hysteresis level is set by THRESH\_LO. This is a level-sensitive gating mode.

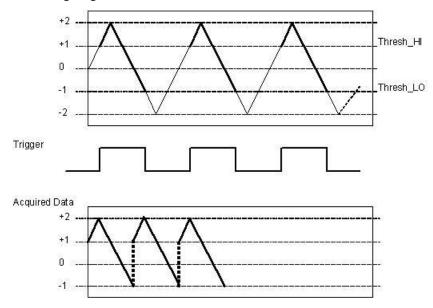


Figure 23. Gate Negative Hysteresis

### **Gate Positive Hysteresis**

Data acquisition is enabled whenever ATRIG goes below the THRESH\_LO level. Acquisition is suspended whenever the ATRIG signal goes above the THRESH\_HI level. The hysteresis level is set by THRESH\_HI. This is a level-sensitive gating mode.

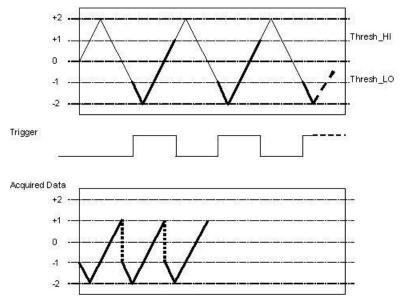


Figure 24. Gate Positive Hysteresis

#### **Gate Inside Window**

Data acquisition is enabled whenever ATRIG is below the THRESH\_HI level and above the THRESH\_LO level. Acquisition is suspended whenever the ATRIG signal is outside of this region. This is a level-sensitive gating mode

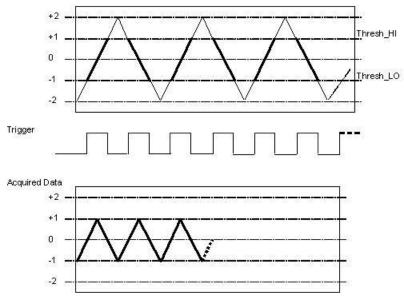


Figure 25. Gate Inside Window

#### **Gate Outside Window**

Data acquisition is enabled whenever ATRIG is above the THRESH\_HI level or below the THRESH\_LO level. Acquisition is suspended whenever the ATRIG signal is between the THRESH\_HI and THRESH\_LO levels. This is a level-sensitive gating mode

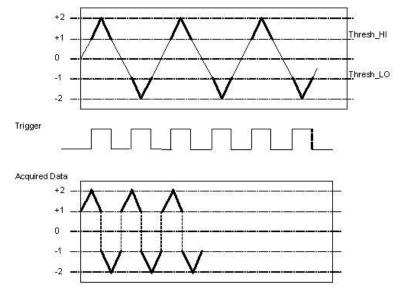


Figure 26. Gate Outside Window

## **Waveform generation timing signals**

The signals that control the timing for the analog output functions on the PCI-DAS6071 are:

- D/A START TRIGGER
- D/A UPDATE
- D/A EXTERNAL TIME BASE

## D/A START TRIGGER signal

The D/A START TRIGGER signal is used to hold off output scans until after a trigger event. The DAQ-Sync "DS D/A START TRIGGER" input or any AUXIN pin can be programmed to serve as the D/A START TRIGGER signal. It is also available as an output on any AUXOUT pin.

When used as an input, the D/A START TRIGGER signal may be software selected as either a positive or negative edge trigger. The selected edge of the D/A START TRIGGER signal causes the DACs to start generating the output waveform.

The D/A START TRIGGER signal can be used as an output to monitor the trigger that initiates waveform generation. The output is an active-high pulse having a width of 50 ns.

Figure 27 and Figure 28 show the input and output timing requirements for the D/A START TRIGGER signal.

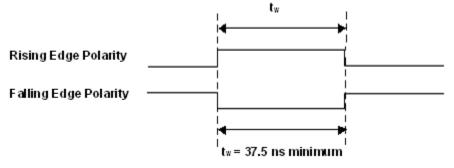


Figure 27. D/A START TRIGGER input signal timing

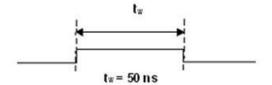


Figure 28. D/A START TRIGGER output signal timing

#### D/A CONVERT signal

The D/A CONVERT signal causes a single output update on the D/A converters. You can program the DAQ-Sync DS D/A UPDATE input or any AUXIN pin to accept the D/A CONVERT signal. It is also available as an output on any AUXOUT pin.

The D/A CONVERT input signal polarity is software selectable. DAC outputs update within 100ns of the selected edge. The D/A CONVERT pulses should be no less than 100  $\mu$ s apart.

When used as an output, the D/A CONVERT signal may be used to monitor the pacing of the output updates. The output has a pulse width of 225 ns with selectable polarity. Figure 29 and Figure 30 show the input and output timing requirements for the D/A CONVERT signal.

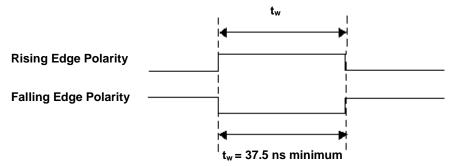


Figure 29. D/A CONVERT input signal timing

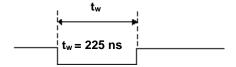


Figure 30. D/A CONVERT output signal timing

## D/A EXTERNAL TIME BASE signal

The D/A EXTERNAL TIME BASE signal can serve as the source for the on-board DAC pacer circuit rather than using the internal time base. Any AUXIN pin can be set programmatically as the source for this signal. The polarity is programmable.

The maximum frequency for the D/A EXTERNAL TIME BASE signal is 20 MHz. The minimum pulse width is 23 ns high or low. There is no minimum frequency specification.

Figure 31 shows the timing requirements for the D/A EXTERNAL TIME BASE signal.

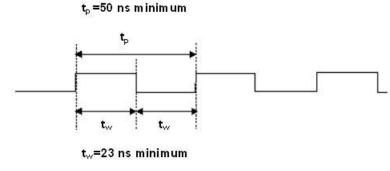


Figure 31. D/A EXTERNAL TIME BASE signal timing

## **General-purpose counter signal timing**

The general-purpose counter signals are:

- CTR1 CLK
- CTR1 GATE
- CTR1 OUT
- CTR2 CLK
- CTR2 GATE
- CTR2 OUT

## CTR1 CLK signal

The CTR1 CLK signal can serve as the clock source for independent user counter 1. It can be selected through software at the CTR1 CLK pin rather than using the on-board 10 MHz or 100 kHz sources. It is also polarity programmable. The maximum input frequency is 10 MHz. There is no minimum frequency specified.

Figure 32 shows the timing requirements for the CTR1 CLK signal.

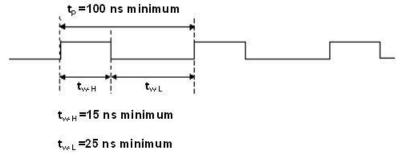


Figure 32. CTR1 CLK signal timing

## CTR1 GATE signal

You can use the CTR1 GATE signal for starting and stopping the counter, saving counter contents, etc. It is polarity programmable and is available at the CTR1 GATE pin.

Figure 33 shows the minimum timing requirements for the CTR1 GATE signal.

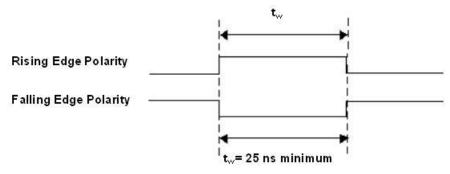


Figure 33. CTR1 GATE signal timing

## CTR1 OUT signal

This signal is present on the CTR1 OUT pin. The CTR1 OUT signal is the output of one of the two user's counters in an industry-standard 82C54 chip. For detailed information on counter operations, refer to the data sheet at <a href="https://www.mccdaq.com/PDFmanuals/82C54.pdf">www.mccdaq.com/PDFmanuals/82C54.pdf</a>.

Figure 34 shows the timing requirements for the CTR1 OUT signal for counter mode 0 and mode 2.

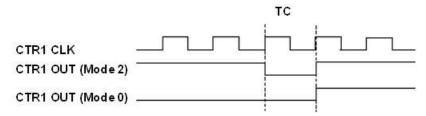


Figure 34. CTR1 OUT signal timing

## CTR2 CLK signal

The CTR2 CLK signal can serve as the clock source for independent user counter 2. It can be selected through software at the CTR2 CLK pin rather than using the on-board 10 MHz or 100 kHz sources. It is also polarity programmable. The maximum input frequency is 10 MHz. There is no minimum frequency specified. Figure 35 shows the timing requirements for the CTR2 CLK signal.

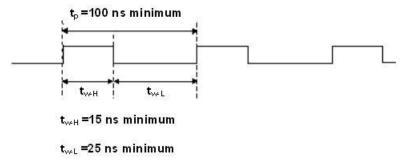


Figure 35. CTR2 CLK signal timing

## **CTR2 GATE signal**

You can use the CTR2 GATE signal for starting and stopping the counter, saving counter contents, etc. It is polarity programmable and is available at the CTR2 GATE pin. Figure 36 shows the timing requirements for the CTR2 GATE signal.

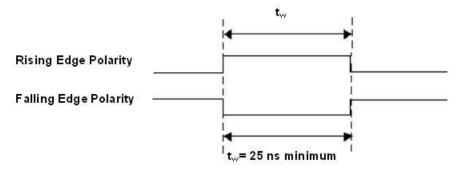


Figure 36. CTR2 GATE signal timing

## **CTR2 OUT signal**

This signal is present on the CTR2 OUT pin. The CTR2 OUT signal is the output of one of the two user's counters in an industry-standard 82C54 chip.

Figure 37 shows the timing of the CTR1 OUT signal for mode 0 and for mode 2.

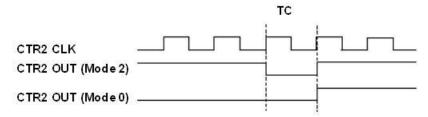


Figure 37. CTR2 OUT signal timing

For detailed information on counter operations, refer to the data sheet at <a href="https://www.mccdaq.com/PDFmanuals/82C54.pdf">www.mccdaq.com/PDFmanuals/82C54.pdf</a>.

## **Calibrating the Board**

## Introduction

You should calibrate the board (using the *InstaCal* utility) after the board has fully warmed up. The recommended warm-up time is 15 minutes. For best results, calibrate the board immediately before making critical measurements. The high resolution analog components on the board are somewhat sensitive to temperature. Pre-measurement calibration ensures that your board is operating with the optimum calibration values.

## Calibration theory

Analog inputs are calibrated for offset and gain. Offset calibration for the analog inputs is performed directly on the programmable gain input amplifier (PGIA), with coarse and fine trim DACs acting on the amplifier.

For input gain calibration, a precision calibration reference is used with coarse and fine trim DACs acting on the ADC (see Figure 38).

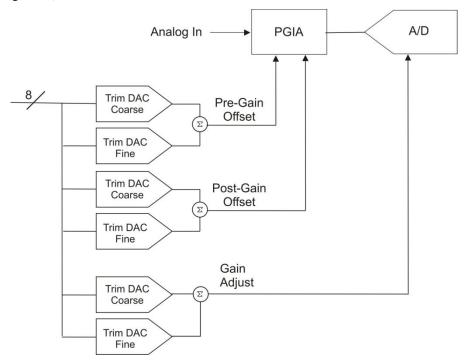


Figure 38. Analog input calibration - basic elements

PCI-DAS6071 User's Guide Calibrating the Board

A similar method is used to calibrate the analog output components. A trim DAC is used to adjust the gain of the DAC. A separate DAC is used to adjust offset on the final output amplifier. The calibration circuits are duplicated for both analog outputs (see Figure 39).

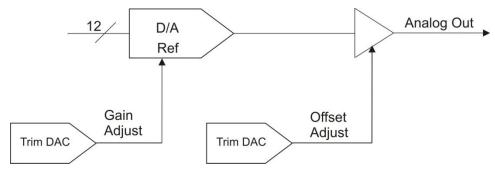


Figure 39. Analog output calibration – basic elements

## **Specifications**

All specifications are subject to change without notice. Typical for 25 °C unless otherwise specified. Specifications in *italic text* are guaranteed by design.

## **Analog input**

Table 1. Analog input specifications

Parameter	Specification				
A/D converter type	Successive approximation				
Maximum sample rate	1.25 MS/s				
Resolution	12 bits, 1-in-4096				
Number of channels	64 single ended / 32 differen	ntial, software selectable			
Input ranges	Bipolar: ±10 V, ±5 V, ±2.5 V, ±1 V, ±0.5 V, ±0.25 V, ±0.1 V, ±0.05 V  Unipolar: 0 to 10 V, 0 to 5 V, 0 to 2 V, 0 to 1 V, 0 to 0.5 V, 0 to 0.2 V,  0 to 0.1 V  Software selectable				
A/D pacing (software programmable)	<ul><li>Internal 40 MHz, 50 ppm</li><li>External Source via AUX</li></ul>	Internal counter – ASIC. Software selectable time base:  Internal 40 MHz, 50 ppm stability  External Source via AUXIN<5:0>, software selectable.			
	External convert strobe: A/I Software paced	O CONVERT			
Burst mode	Software selectable option,	burst rate = $800 \text{ nS}$			
A/D gate sources	External digital:	A/D GATE			
	External analog:	ATRIG input CH0 IN through CH63 IN			
A/D gating modes	External digital:	Programmable, active high or active low, level or edge			
	External analog:	Refer to Analog trigger on page 41			
A/D trigger sources	External digital:	A/D START TRIGGER A/D STOP TRIGGER			
	External analog:	ATRIG input CH0 IN through CH63 IN			
A/D triggering modes	External digital:	Software-configurable for rising or falling edge.			
	External analog:	Refer to Analog trigger on page 41			
	Pre-/Post-trigger:	Unlimited number of pre-trigger samples, 16 Meg post-trigger samples.			
ADC pacer out	Available at user connector:	: A/D PACER OUT			
RAM buffer size	8K samples				
Data transfer	DMA Programmed I/O				
DMA modes		Demand or non-demand using scatter gather.			
Configuration memory	Up to 8 K elements. Programmable channel, gain, and offset				
Streaming-to-disk rate	1.25 MS/s, system dependent				

PCI-DAS6071 User's Guide Specifications

## **Accuracy**

1.25 MS/s rate, single channel operation and a 15-minute warm-up. Accuracies listed are for measurements made following an internal calibration. They are valid for operational temperatures within  $\pm 1^{\circ}$  C of internal calibration temperature and  $\pm 10^{\circ}$ C of factory calibration temperature. Calibrator test source high side tied to channel 0 high and low side tied to channel 0 low. Low-level ground is tied to channel 0 low at the user connector.

Table 2. Absolute accuracy specifications

Range	Absolute Accuracy (LSB)
±10V	±2.9
±5V	±2.1
±2.5V	±3.0
±1V	±3.0
±500mV	±3.0
±250mV	±3.1
±100mV	±3.3
±50mV	±3.7
0 to 10V	±2.8
0 to 5V	±4.4
0 to 2V	±4.4
0 to 1V	±4.5
0 to 500mV	±4.6
0 to 200mV	±4.8
0 to 100mV	±5.2

Table 3. Absolute accuracy components specifications – all values are (±)

Range	% of	Offset	Noise + Qua	ntization (mV)	Temp Drift (%/DegC)	Absolute
	Reading	(mV)	Single Pt	Averaged (Note 1)		Accuracy at FS (mV)
±10V	0.0714	6.38	6.10	0.846	0.0010	14.369
±5V	0.0314	3.20	3.05	0.423	0.0005	5.193
±2.5V	0.0714	1.61	1.53	0.211	0.0010	3.605
±1V	0.0714	0.653	0.610	0.085	0.0010	1.452
±500mV	0.0714	0.335	0.305	0.042	0.0010	0.735
±250mV	0.0714	0.176	0.208	0.024	0.0010	0.379
±100mV	0.0714	0.081	0.098	0.011	0.0010	0.163
±50mV	0.0714	0.049	0.071	0.007	0.0010	0.091
0 to 10V	0.0314	3.20	3.05	0.423	0.0005	6.765
0 to 5V	0.0714	1.61	1.53	0.211	0.0010	5.391
0 to 2V	0.0714	0.653	0.610	0.085	0.0010	2.167
0 to 1V	0.0714	0.335	0.305	0.042	0.0010	1.092
0 to 500mV	0.0714	0.176	0.208	0.024	0.0010	0.558
0 to 200mV	0.0714	0.081	0.098	0.011	0.0010	0.235
0 to 100mV	0.0714	0.049	0.071	0.007	0.0010	0.127

Note 1: Averaged measurements assume dithering and averaging of 100 single-channel readings.

Each PCI-DAS6071 is tested at the factory to assure the board's overall error does not exceed absolute accuracy limits described in Table 2.

PCI-DAS6071 User's Guide Specifications

Table 4. Relative accuracy specifications - all values are (±)

Range	Relative Accuracy (mV)		
	Single Point	Averaged (Note 2)	
±10V	7.37	1.11	
±5V	3.68	0.557	
±2.5V	1.84	0.278	
±1V	0.737	0.111	
±500mV	0.368	0.056	
±250mV	0.238	0.032	
±100mV	0.111	0.015	
±50mV	0.082	0.009	
0 to 10V	3.68	0.557	
0 to 5V	1.84	0.278	
0 to 2V	0.737	0.111	
0 to 1V	0.368	0.056	
0 to 500mV	0.238	0.032	
0 to 200mV	0.111	0.015	
0 to 100mV	0.082	0.009	

**Note 2:** Averaged measurements assume dithering and averaging of 100 single-channel readings.

Relative accuracy is defined as the measured deviation from a straight line drawn between measured endpoints of the transfer function. ADC resolution, noise and front-end non-linearity are included in this measurement.

Table 5. Differential non-linearity specifications

All ranges	+0.5 I CD tvm	+1 O I CD may
All ranges	±0.5 LSB typ	±1.0 LSB max

### **Settling Time**

Settling time is defined here as the time required for a channel to settle to within a specified accuracy in response to a full-scale (FS) step. Two channels are scanned at a specified rate. A -FS DC signal is presented to channel 1; a +FS DC signal is presented to channel 0.

Table 6. Settling time specifications

Condition	Range	Accuracy			
			±0.012% (±0.5 LSB)	±0.024% (±1.0 LSB)	±0.098% (±4.0 LSB)
Same range to same range	±10 V	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.9 μS 2.0 μS
sume range	±5 V	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.5 μS
	±2.5 V	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.3 μS
	±1 V	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.3 μS
	±500 mV	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.3 μS
	±250 mV	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.3 μS
	±100 mV	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.3 μS
	±50 mV	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.5 μS
	0 to 10 V	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.5 μS
	0 to 5 V	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.3 μS
	0 to 2 V	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.3 μS
	0 to 1 V	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.3 μS
	0 to 500 mV	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.3 μS
	0 to 200 mV	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.3 μS
	0 to 100 mV	Typ Max	3.0 μS 5.0 μS	1.9 μS 2.5 μS	1.2 μS 1.5 μS

#### **Parametrics**

Table 7. Parametric specifications

Parameter	Specification	
Max working voltage (signal + common-mode)	Input must remain within ±11 V of ground	
CMRR @ 60 Hz	±10 V: 95 dB ±5 V, 0 to 10 V: 100 dB All other ranges: 106 dB	
Small signal bandwidth, all ranges	1.6 MHz	
Large signal bandwidth, all ranges	1.0 MHz	
Input coupling	DC	
Input impedance	$100~G\Omega$ in parallel with $100pF$ in normal operation.	
Input bias current	±200 pA	
Input offset current	±100 pA	
Absolute maximum input voltage	Power ON: ±25 V, Power OFF: ±15 V (±20 mA Note 3) Protected inputs:  • CH<63:0> IN  • AISENSE	
Power on and reset state	CH0 IN, single-ended mode, 0 V to 0.1 V input range (Note 4)	
Crosstalk, DC to 100 kHz	Adjacent channels: -75 dB	
	All other channels: –90 dB	

**Note 3:** The analog input sink/source current must be limited to an maximum of  $\pm 20$  mA in the power OFF state to prevent damage to the board. A 1000  $\Omega$  (½ W) current limiting resistor should be placed in series with each analog input channel being used in applications where the power OFF state sink/source current into the board can exceed  $\pm 20$  mA. Resistance values >1000  $\Omega$  may adversely affect the noise and settling time performance of the board.

**Note 4:** Care should be taken to avoid the application of an input voltage to CH0 IN that could overdrive the analog input circuit. Any unused analog input channel should be connected to LLGND.

#### **Noise Performance**

Table 8 summarizes the noise performance for the PCI-DAS6071. Noise distribution is determined by gathering 50K samples with inputs tied to ground at the user connector. Samples are gathered at the maximum specified single-channel sampling rate. Specification applies to both single-ended and differential modes of operation.

Table 8. Analog input noise performance specifications (not including quantization)

Range	Counts Dithered	LSBrms Dithered	Counts Undithered	LSBrms Undithered
±10 V	5	0.5	3	0.25
±5 V	5	0.5	3	0.25
±2.5 V	5	0.5	3	0.25
±1 V	5	0.5	3	0.25
±500 mV	6	0.5	3	0.25
±250 mV	6	0.6	4	0.4
±100 mV	7	0.7	5	0.5
±50 mV	9	0.9	8	0.8
0 to 10 V	5	0.5	3	0.25
0 to 5 V	5	0.5	3	0.25
0 to 2 V	5	0.5	3	0.25
0 to 1 V	6	0.5	3	0.25
0 to 500 mV	6	0.6	4	0.4
0 to 200 mV	7	0.7	5	0.5
0 to 100 mV	9	0.9	8	0.8

### **Analog output**

Table 9. AO specifications

Parameter	Specification
D/A converter type	Double-buffered, multiplying
Resolution	12-bits, 1-in-4096
Number of channels	2 voltage output
Voltage range	±10 V, 0 to 10 V, ±EXT REF., 0 to EXT REF., software selectable
Monotonicity	12-bits, guaranteed
Slew rate	20 V/μs min
Settling time (full scale step)	$3.0 \mu S$ to $\pm 0.5 LSB$ accuracy
Noise	200 μVrms, DC to 1MHz BW
Glitch energy	±20 mV @ 1.5 μS duration measured at mid-scale transition.
Current drive	±5 mA
Output short-circuit duration	Indefinite @ 25 mA
Output coupling	DC
Output impedance	$0.1~\Omega$ max
Gain temperature coefficient, internal or external reference	25 ppm/°C
Offset temperature coefficient	±50 μV/°C
Power up and reset	DACs cleared to 0 volts ±200 mV max

Table 10. Analog output absolute accuracy specifications

Range	Absolute Accuracy (LSB)
±10 V	±1.7 LSB
0 to 10 V	±2.3 LSB

Table 11. Analog output absolute accuracy components specifications

Range	% of Reading	Offset (mV)	Temp Drift (%/DegC)	Absolute Accuracy at FS (mV)
±10 V	±0.0219	±5.93	±0.0005	±8.127
0 to 10 V	±0.0219	±3.49	±0.0005	±5.685

Each PCI-DAS6071 is tested at the factory to assure the board's overall error does not exceed the limits listed in Table 10.

Table 12. Relative accuracy specifications

Range	Relative Accuracy	
All ranges	±0.3 LSB, typ	±0.5 LSB, max

Relative accuracy is defined as the measured deviation from a straight line drawn between measured endpoints of the transfer function.

Table 13. Differential non-linearity specifications

|--|

### Analog output pacing and triggering

Table 14. External reference input (D/A EXTREF)specifications

Parameter	Specification		
DAC pacing	Internal counter – ASIC. Selectable time base:		
(software programmable)	■ Internal 40 MHz		
	<ul> <li>External Source via AUX</li> </ul>	IN<5:0>, SW selectable.	
	External convert strobe:	D/A UPDATE	
	Software paced		
DAC gate source	External digital:	D/A START TRIGGER	
(software programmable)	External analog:	ATRIG input	
		CH0 IN through CH63 IN	
	Software gated		
DAC gating modes	External digital:		
	<ul> <li>Programmable, active high or active low, level or edge</li> </ul>		
	External analog:	Refer to Analog trigger on page 41	
DAC trigger sources	External digital:	D/A START TRIGGER	
	External analog:	ATRIG input	
		CH0 IN through CH63 IN	
	Software triggered		
DAC triggering modes	External digital:		
	<ul> <li>Software-configurable for rising or falling edge.</li> </ul>		
	External analog:	Refer to Analog trigger on page 41	
DAC pacer out	Available at user connector:	D/A PACER OUT	
RAM buffer size	16 K samples		
Data transfer	DMA		
	Programmed I/O		
	Update DACs individually or simultaneously, software selectable.		
DMA modes	Demand or Non-Demand using scatter gather.		
Waveform generation throughput	1 MS/s max per channel, 2 c	channels simultaneous	

### Analog output external reference input (D/A EXTREF)

Table 15. External reference input (D/A EXTREF) specifications

Parameter	Specification
Range	±11 V
Overvoltage protection	±25 V powered on, ±15V powered off
Input impedance	10 kΩ
Bandwidth (-3 dB)	1 MHz
Gain Error – EXTREF mode	0 to 0.5%; not adjustable.

# **Analog trigger**

Table 16. Analog trigger specifications

Parameter	Specification	
Analog trigger sources	External: ATRIG input	
Software selectable	CH0 IN through CH63 IN, first channel in scan	
Analog trigger levels	ATRIG input: ±10 V	
	CH0 IN through CH63 IN: ± Full-scale, range dependent	
Analog trigger modes	External analog: Software-configurable for:	
	Positive or Negative slope	
Analog gate modes	External analog: Software-configurable for:	
	Above or below reference	
	<ul> <li>Positive or negative hysteresis</li> </ul>	
	■ In or out of window	
Resolution	8-bits, 1-in-256	
Accuracy	±5% Full-scale range max	
Bandwidth (–3 dB)	ATRIG input: 1.3 MHz	
	CH0 IN through CH63 IN: 2.0 MHz	

# Analog I/O calibration

Table 17. Analog I/O calibration specifications

Parameter	Specification	
Recommended warm-up time	15 minutes	
Calibration	Auto-calibration, calibration factors for each range stored on board in non-volatile RAM.	
Onboard calibration reference	DC Level: 5.000 V± 2.5 mv. Actual measured values stored in EEPROM.	
	Tempco: 5ppm/°C max, 2 ppm/°C typ	
	Long-term stability: 20 ppm, T = 1000 hrs, non-cumulative	
Calibration interval	1 year	

# Digital I/O

Table 18. Digital I/O specifications

Parameter	Specification	
Digital type	Discrete, 5V/TTL compatible	
Number of I/O	8	
Configuration	8 bits, independently programmable for input or output. All pins pulled up to +5 V via 47 K resistors (default). Positions available for pull-down to ground. Hardware selectable via solder gap.	
Input high voltage	2.0 V min, 7.0 V absolute max	
Input low voltage	0.8 V max, -0.5 V absolute min	
Output high voltage (IOH = -32 mA)	3.80 V min, 4.20 V typ	
Output low voltage (IOL = 32 mA)	0.55 V max, 0.22 V typ	
Data transfer	Programmed I/O	
Power-up / reset state	Input mode (high impedance)	

# Interrupts

Table 19. Interrupt specifications

Parameter	Specification		
Interrupts	PCI INTA# - mapped to IRQn via PCI BIOS at boot-time		
Interrupt enable	Programmable through PLX9080		
ADC interrupt sources	DAQ_ACTIVE: Interrupt is generated when a DAQ sequence is active.		
(software programmable)	DAQ_STOP:	Interrupt is generated when A/D Stop Trigger In is detected.	
	DAQ_DONE:	Interrupt is generated when a DAQ sequence completes.	
	DAQ_FIFO_1/4_	=	
	D. C. CD. CT. D.	Interrupt is generated when ADC FIFO is ¼ full.	
	DAQ_SINGLE:	Interrupt is generated after each conversion completes.	
	DAQ_EOSCAN: Interrupt is generated after the last channel is converted in multi-channel scans.		
	DAQ_EOSEQ:	Interrupt is generated after each interval delay during multi- channel scans.	
DAC interrupt sources (software programmable)	DAC_ACTIVE:	Interrupt is generated when DAC waveform circuitry is active.	
	DAC_DONE:	Interrupt is generated when a DAC sequence completes.	
	DAC_FIFO_1/4_EMPTY:		
		Interrupt is generated DAC FIFO is ¼ empty.	
	DAC_HIGH_CH	ANNEL:	
		Interrupt is generated when the DAC high channel output is	
		updated.	

### **Counters**

Table 20. Counter specifications

Parameter	Specification	
User counter type	82C54	
Number of channels	2	
Resolution	16-bits	
Compatibility	5V/TTL	
CTRn base clock source (software selectable)	Internal 10 MHz, Internal 100 kHz, or External connector (CTRn CLK)	
Internal 10 MHz clock source stability	50 ppm	
Counter n gate	Available at connector (CTRn GATE)	
Counter n output	Available at connector (CTRn OUT)	
Clock input frequency	10 MHz max	
High pulse width (clock input)	15 ns min	
Low pulse width (clock input)	25 ns min	
Gate width high	25 ns min	
Gate width low	25 ns min	
Input low voltage	0.8 V max	
Input high voltage	2.0 V min	
Output low voltage	0.4 V max	
Output high voltage	3.0 V min	

### Configurable AUXIN<5:0>, AUXOUT<2:0> external trigger/clocks

The PCI-DAS6071 provides nine user-configurable trigger/clock pins available at the 100-pin I/O connector. Of these, six are configurable as inputs while three are configurable as outputs.

Table 21. Configurable triggers/clocks specifications

Parameter	Specification		
AUXIN<5:0> sources (software	A/D CONVERT:	External ADC convert strobe	
selectable)	A/D TIMEBASE IN:	External ADC pacer timebase	
	A/D START TRIGGER:	ADC Start Trigger	
	A/D STOP TRIGGER:	ADC Stop Trigger	
	A/D PACER GATE:	External ADC gate	
	D/A START TRIGGER:	DAC trigger/gate	
	D/A UPDATE:	DAC update strobe	
	D/A TIMEBASE IN:	External DAC pacer time base	
AUXOUT<2:0> sources (software	STARTSCAN:	A pulse indicating start of conversion	
selectable)	SSH:	Active signal that terminates at the start of the last	
		conversion in a scan	
	A/D STOP:	Indicates end of scan	
	A/D CONVERT:	ADC convert pulse	
	SCANCLK:	Delayed version of ADC convert	
	CTR1 CLK:	CTR1 clock source	
	D/A UPDATE:	D/A update pulse	
	CTR2 CLK:	CTR2 clock source	
	A/D START TRIGGER:	ADC Start Trigger Out	
	A/D STOP TRIGGER:	ADC Stop Trigger Out	
	A/D PACER GATE:	External ADC gate	
	D/A START TRIGGER:	DAC Start Trigger Out	
Default selections:	AUXIN0:	A/D CONVERT	
	AUXIN1:	A/D START TRIGGER	
	AUXIN2:	A/D STOP TRIGGER	
	AUXIN3:	D/A UPDATE	
	AUXIN4:	D/A START TRIGGER	
	AUXIN5:	A/D PACER GATE	
	AUXOUT0:	D/A UPDATE	
	AUXOUT1:	A/D CONVERT	
	AUXOUT2:	SCANCLK	
Compatibility	5 V/TTL		
Edge-sensitive polarity	Rising/falling, software selectable		
Level-sensitive polarity	Active high/active low, software selectable		
Minimum pulse width	37.5 nS		

## **DAQ-Sync inter-board triggers/clocks**

The DAQ-Sync bus provides inter-board triggering and synchronization capability. Five trigger/strobe I/O pins and one clock I/O pin are provided on a 14-pin header. The DAQ-Sync signals use dedicated pins. Only the direction may be set.

Table 22. DAQ-Sync signal specifications

Connector	Signal name
DAQ-Sync	DS A/D START TRIGGER
	DS A/D STOP TRIGGER
	DS A/D CONVERT
	DS D/A UPDATE
	DS D/A START TRIGGER
	SYNC CLK

### **Power consumption**

Table 23. Power consumption specifications

Parameter	Specification	
+5 V	0.9 A typ, 1.1 A max. Does not include power consumed through the I/O connector	
+5 V available at I/O connector	r 1 A max, protected with a resettable fuse	

#### **Environmental**

Table 24. Environmental specifications

Parameter	Specification	
Operating temperature range	0 °C to 55 °C	
Storage temperature range		
Humidity	0% to 90% non-condensing	

#### Mechanical

Table 25. Mechanical specifications

Parameter	Specification	
Card dimensions (L × W × H) PCI half card: 174.4 (6.87) × 106.9 (4.21) × 11.65 mm (0.46 in.)		

# **DAQ-Sync connector**

Table 26. DAQ-Sync connector specifications

Parameter	Specification	
Connector type	14-pin right-angle 100 mil box header	
Compatible cables	MCC p/n: CDS-14-x, 14-pin ribbon cable. x = number of boards (2 to 5)	

Table 27. DAQ-Sync connector pinout

Pin	Signal Name
1	DS A/D START TRIGGER
2	GND
3	DS A/D STOP TRIGGER
4	GND
5	DS A/D CONVERT
6	GND
7	DS D/A UPDATE
8	GND
9	DS D/A START TRIGGER
10	GND
11	RESERVED
12	GND
13	SYNC CLK
14	GND

#### **SCSI** connector

Table 28. SCSI connector specifications

Parameter	Specification	
Connector type	Shielded SCSI 100 D-type	
Compatible Cables	C100HD50-x, unshielded ribbon cable. $x = 3$ or 6 feet	
	C100MMS-x, shielded round cable.	
	x = 1, 2,  or  3  meters	
Compatible accessory products   ISO-RACK16/P		
(with the C100HD50-x cable)	ISO-DA02/P	
	BNC-16SE	
	BNC-16DI	
	CIO-MINI50	
	CIO-TERM100	
	SCB-50	
Compatible accessory products	SCB-100	
(with the C100MMS-x cable)		

Table 29. 32-channel differential mode pinout

Pin	Signal Name	Pin	Signal Name
1	LLGND	51	LLGND
2	CH0 IN HI	52	CH16 IN HI
3	CH0 IN LO	53	CH16 IN LO
4	CH1 IN HI	54	CH17 IN HI
5	CH1 IN LO	55	CH17 IN LO
6	CH2 IN HI	56	CH18 IN HI
7	CH2 IN LO	57	CH18 IN LO
8	CH3 IN HI	58	CH19 IN HI
9	CH3 IN LO	59	CH19 IN LO
10	CH4 IN HI	60	CH20 IN HI
11	CH4 IN LO	61	CH20 IN LO
12	CH5 IN HI	62	CH21 IN HI
13	CH5 IN LO	63	CH21 IN LO
14	CH6 IN HI	64	CH22 IN HI
15	CH6 IN LO	65	CH22 IN LO
16	CH7 IN HI	66	CH23 IN HI
17	CH7 IN LO	67	CH23 IN LO
18	LLGND	68	LLGND
19	CH8 IN HI	69	CH24 IN HI
20	CH8 IN LO	70	CH24 IN LO
21	CH9 IN HI	71	CH25 IN HI
22	CH9 IN LO	72	CH25 IN LO
23	CH10 IN HI	73	CH26 IN HI
24	CH10 IN LO	74	CH26 IN LO
25	CH11 IN HI	75	CH27 IN HI
26	CH11 IN LO	76	CH27 IN LO
27	CH12 IN HI	77	CH28 IN HI
28	CH12 IN LO	78	CH28 IN LO
29	CH13 IN HI	79	CH29 IN HI
30	CH13 IN LO	80	CH29 IN LO
31	CH14 IN HI	81	CH30 IN HI
32	CH14 IN LO	82	CH30 IN LO
33	CH15 IN HI	83	CH31 IN HI
34	CH15 IN LO	84	CH31 IN LO
35	AISENSE	85	DIO0
36	D/A OUT 0	86	DIO1
37	D/A GND	87	DIO2
38	D/A OUT1	88	DIO3
39	PC +5 V	89	DIO4
40	AUXOUT0 / D/A PACER OUT	90	DIO5
41	AUXOUT1 / A/D PACER OUT	91	DIO6
42	AUXOUT2 / SCANCLK	92	DIO7
43	AUXINO / A/D CONVERT / ATRIG	93	CTR1 CLK
44	D/A EXTREF	94	CTR1 GATE
45	AUXIN1 / A/D START TRIGGER	95	CTR1 OUT
46	AUXIN2 / A/D STOP TRIGGER	96	GND
47	AUXIN3 / D/A UPDATE	97	CTR2 CLK
48	AUXIN4 / D/A START TRIGGER	98	CTR2 GATE
49	AUXIN5 / A/D PACER GATE	99	CTR2 OUT
50	GND	100	GND
	J.10	100	J., D

Table 30. 64-channel single-ended mode pinout

Pin	Signal Name	Pin	Signal Name
1	LLGND	51	LLGND
2	CH0 IN	52	CH16 IN
3	CH32 IN	53	CH48 IN
4	CH1 IN	54	CH17 IN
5	CH33 IN	55	CH49 IN
6	CH2 IN	56	CH18 IN
7	CH34 IN	57	CH50 IN
8	CH3 IN	58	CH19 IN
9	CH35 IN	59	CH51 IN
10	CH4 IN	60	CH20 IN
11	CH36 IN	61	CH52 IN
12	CH5 IN	62	CH21 IN
13	CH37 IN	63	CH53 IN
14	CH6 IN	64	CH22 IN
15	CH38 IN	65	CH54 IN
16	CH7 IN	66	CH23 IN
17	CH39 IN	67	CH55 IN
18	LLGND	68	LLGND
19	CH8 IN	69	CH24 IN
20	CH40 IN	70	CH56 IN
21	CH9 IN	71	CH25 IN
22	CH41 IN	72	CH57 IN
23	CH10 IN	73	CH26 IN
24	CH42 IN	74	CH58 IN
25	CH11 IN	75	CH27 IN
26	CH43 IN	76	CH59 IN
27	CH12 IN	77	CH28 IN
28	CH44 IN	78	CH60 IN
29	CH13 IN	79	CH29 IN
30	CH45 IN	80	CH61 IN
31	CH14 IN	81	CH30 IN
32	CH46 IN	82	CH62 IN
33	CH15 IN	83	CH31 IN
34	CH47 IN	84	CH63 IN
35	AISENSE	85	DIO0
36	D/A OUT 0	86	DIO1
37	D/A GND	87	DIO2
38	D/A OUT1	88	DIO3
39	PC +5 V	89	DIO4
40	AUXOUT0 / D/A PACER OUT	90	DIO5
41	AUXOUT1 / A/D PACER OUT	91	DIO6
42	AUXOUT2 / SCANCLK	92	DIO7
43	AUXIN0 / A/D CONVERT / ATRIG	93	CTR1 CLK
44	D/A EXTREF	94	CTR1 GATE
45	AUXIN1 / A/D START TRIGGER	95	CTR1 OUT
46	AUXIN2 / A/D STOP TRIGGER	96	GND
47	AUXIN3 / D/A UPDATE	97	CTR2 CLK
48	AUXIN4 / D/A START TRIGGER	98	CTR2 GATE
49	AUXIN5 / A/D PACER GATE	99	CTR2 OUT
50	GND	100	GND
	1 - 7	1	<u> </u>

# CE Declaration of Conformity

Manufacturer: Measurement Computing Corporation

Address: 10 Commerce Way

Suite 1008

Norton, MA 02766

USA

Category: Electrical equipment for measurement, control and laboratory use.

Measurement Computing Corporation declares under sole responsibility that the product

#### PCI-DAS6071

to which this declaration relates is in conformity with the relevant provisions of the following standards or other documents:

EC EMC Directive 2004/108/EC: General Requirements, EN 61326-1:2006 (IEC 61326-1:2005).

#### **Emissions:**

- EN 55011 (2007) / CISPR 11(2003): Radiated emissions: Group 1, Class A
- EN 55011 (2007) / CISPR 11(2003): Conducted emissions: Group 1, Class A

Immunity: EN 61326-1:2006, Table 3.

- IEC 61000-4-2 (2001): Electrostatic Discharge immunity.
- IEC 61000-4-3 (2002): Radiated Electromagnetic Field immunity.
- IEC 61000-4-4 (2004): Electric Fast Transient Burst Immunity.
- IEC 61000-4-5 (2001): Surge Immunity.
- IEC 61000-4-6 (2003): Radio Frequency Common Mode Immunity.
- IEC 61000-4-11 (2004): Voltage Interrupts.

To maintain compliance to the standards of this declaration, the following conditions must be met.

- The host computer, peripheral equipment, power sources, and expansion hardware must be CE compliant.
- All I/O cables must be shielded, with the shields connected to ground.
- I/O cables must be less than 3 meters (9.75 feet) in length.
- The host computer must be properly grounded.
- Equipment must be operated in a controlled electromagnetic environment as defined by Standards EN 61326-1:2006, or IEC 61326-1:2005.

Declaration of Conformity based on tests conducted by Chomerics Test Services, Woburn, MA 01801, USA in May, 2004. Test records are outlined in Chomerics Test Report #EMI3889.04. Further testing was conducted by Chomerics Test Services, Woburn, MA. 01801, USA in January, 2009. Test records are outlined in Chomerics Test report #EMI5243.09.

We hereby declare that the equipment specified conforms to the above Directives and Standards.

Carl Haapaoja, Director of Quality Assurance

Call taggage

Measurement Computing Corporation 10 Commerce Way

**Suite 1008** 

Norton, Massachusetts 02766

(508) 946-5100

Fax: (508) 946-9500

E-mail: info@mccdaq.com www.mccdaq.com